



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
4-29-2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/25/14 Cal 1 0.099 AC Cal 2 0.097 AC

Between the dates of 3/28/14 and 4/25/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


4/28/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 28th, day of April, 2014.

Elma Escarcega

Notary Public

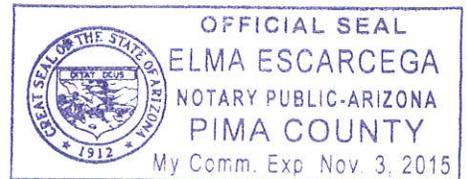

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7 sc
TSY
4-29-2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 4/25/14 TIME 948
INTOXILYZER SERIAL # 80-000331 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
04/25/2014

09:50:30

2 of 7
sc
TS4

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 09:51:32 |
| Diagnostic Test | Pass | 09:52:06 |
| Air Blank | 0.000 | 09:52:33 |
| 0.100 Cal Check | 0.099 | 09:52:54 |
| Air Blank | 0.000 | 09:53:24 |
| Subject Test | 0.000 | 09:53:51 |
| Air Blank | 0.000 | 09:54:19 |
| Wait | ABT* | |
| Air Blank | 0.000 | 09:54:56 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
04/25/2014

09:56:33

307
SC
TS4

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 09:57:30 |
| Diagnostic Test | Pass | 09:58:04 |
| Air Blank | 0.000 | 09:58:31 |
| 0.100 Cal Check | RFI* | 09:58:44 |
| Air Blank | IPS** | 09:58:52 |

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
04/25/2014

09:58:59

4 of 7
sc
B4

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 09:59:54 |
| Diagnostic Test | Pass | 10:00:28 |
| Air Blank | 0.000 | 10:00:56 |
| 0.100 Cal Check | 0.098 | 10:01:17 |
| Air Blank | 0.000 | 10:01:46 |
| Subject Test | 0.369* | 10:04:54 |
| Air Blank | 0.000 | 10:05:27 |
| Air Blank | 0.000 | 10:05:55 |
| Subject Test | XXX** | 10:06:19 |
| Air Blank | 0.000 | 10:06:52 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
04/25/2014 10:09:09

5 of 7
5c
TS4

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:09:37 |
| 0.100 Cal Check | 0.097 | 10:09:58 |
| Air Blank | 0.000 | 10:10:28 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
sc
TBH
4.29.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 4/25/14 TIME 948
Intoxilyzer Serial # 80-000331 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
04/25/2014 09:48:47

7 of 7
S
B4

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 09:49:15 |
| 0.100 Cal Check | 0.099 | 09:49:36 |
| Air Blank | 0.000 | 09:50:06 |

Preventative Maintenance
Performed -- Timer Reset