



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS6
9-29-2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/12/14 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 8/14/14 and 9/12/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


9/15/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 15th day of September, 2014.

Edma Escarcega

Notary Public

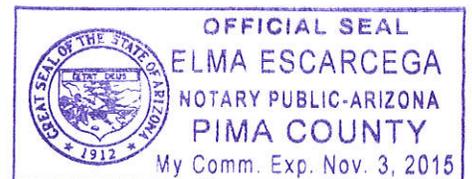

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1077<sub>52</sub>  
TS4  
9.29.2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/12/14 TIME 948

INTOXILYZER SERIAL # 80-000331 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

( ✓ ) 1. Clock time check.

( ✓ ) 2. Date check.

OPERATIONAL TESTS

( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.

( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 095 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
09/12/2014

09:47:05

2 of 7<sub>sc</sub>  
BY

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:48:03
Diagnostic Test	Pass	09:48:37
Air Blank	0.000	09:49:05
0.100 Cal Check	0.096	09:49:26
Air Blank	0.000	09:49:55
Subject Test	0.000	09:50:09
Air Blank	0.000	09:50:36
Wait	ABT*	
Air Blank	0.000	09:51:17

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
09/12/2014

09:51:32

3 of 7  
124

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:52:29
Diagnostic Test	Pass	09:53:03
Air Blank	0.000	09:53:30
0.100 Cal Check	RFI*	09:53:47
Air Blank	IPS**	09:54:15

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
09/12/2014

09:54:33

4 of 7 sc  
Bh

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP

DOB: 11/11/1980

Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:55:27
Diagnostic Test	Pass	09:56:01
Air Blank	0.000	09:56:29
0.100 Cal Check	0.096	09:56:50
Air Blank	0.000	09:57:20
Subject Test	0.310*	10:00:27
Air Blank	0.000	10:01:01
Air Blank	0.000	10:01:28
Subject Test	XXX**	10:01:48
Air Blank	0.000	10:02:22

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
09/12/2014 10:02:38

5 of 7  
sc  
TS4

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:03:07
0.100 Cal Check	0.095	10:03:27
Air Blank	0.000	10:03:57

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
JS  
9.29.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 9/12/14 TIME 948  
Intoxilyzer Serial # 80-000531 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
09/12/2014 09:48:15

7 of 75c  
134

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:48:44
0.100 Cal Check	0.096	09:49:04
Air Blank	0.000	09:49:34

Preventative Maintenance  
Performed -- Timer Reset