



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

Old
 10/23/15

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10.8.2015 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 9.10.2015 and 10.8.2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
 The instrument was returned to service on .

Terry Gallegos

Terry S. Gallegos

10.16.2015


ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 16th, day of October, 2015.

Edma Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

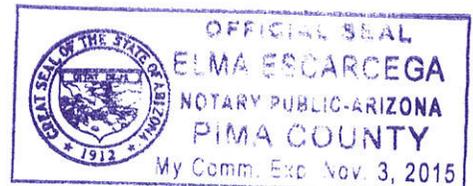

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10 of 7
TS4
Ad
10/23/15

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 10-8-2015 TIME 1201

INTOXILYZER SERIAL # 80-000331 LOCATION DUI-HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
10/08/2015

12:10:05

2 of 7
TS
dd

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:10:52
Diagnostic Test	Pass	12:11:26
Air Blank	0.000	12:11:53
0.100 Cal Check	0.098	12:12:14
Air Blank	0.000	12:12:44
Subject Test	0.000	12:13:00
Air Blank	0.000	12:13:28
Wait	ABT*	
Air Blank	0.000	12:14:02

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
10/08/2015

12:14:41

3 of 7
TSY
JLB

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:15:40
Diagnostic Test	Pass	12:16:14
Air Blank	0.000	12:16:42
0.100 Cal Check	RFI*	12:16:58
Air Blank	IPS**	12:17:09

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 23-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
10/08/2015

12:18:17

4 of 7
TSy
dd

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:19:04
Diagnostic Test	Pass	12:19:38
Air Blank	0.000	12:20:05
0.100 Cal Check	0.098	12:20:26
Air Blank	0.000	12:20:56
Subject Test	0.107*	12:24:03
Air Blank	0.000	12:24:35
Air Blank	0.000	12:25:03
Subject Test	XXX**	12:25:32
Air Blank	0.000	12:26:07

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 23-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
10/08/2015

12:26:20

50f7
TSy
JL

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:26:49
0.100 Cal Check	0.097	12:27:10
Air Blank	0.000	12:27:40

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6077
TS4
all
10/23/15

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 10-8-2015 TIME 1201
Intoxilyzer Serial # 80-000331 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
10/08/2015

12:01:59

7 of 7
TSg
dlb

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:02:28
0.100 Cal Check	0.097	12:02:49
Air Blank	0.000	12:03:19

Preventative Maintenance
Performed -- Timer Reset