



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

*Handwritten:* 11-25-15

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11-5-2015 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 10-8-2015 and 11-5-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

*Handwritten signature:* Terry S. Gallegos

Terry S. Gallegos


11-25-2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 25<sup>th</sup> day of November, 2015.

*Handwritten signature:* Elma Escarcega

Notary Public

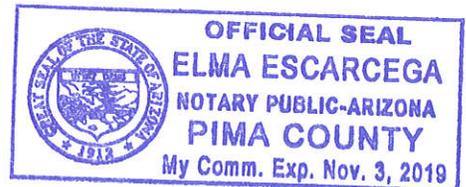

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10F8  
TS4  
JWA  
11-25-15

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 11-5-2015 TIME 1550

INTOXILYZER SERIAL # 80-000331 LOCATION DUI-HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
11/05/2015

15:54:05

2018  
Dh  
JL

---

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	15:54:57
Diagnostic Test	Pass	15:55:31
Air Blank	0.000	15:55:59
0.100 Cal Check	0.099	15:56:19
Air Blank	0.000	15:56:49
Subject Test	0.000	15:57:15
Air Blank	0.000	15:57:43
Wait	ABT*	
Air Blank	0.000	15:58:16

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
11/05/2015 16:01:10

30f8  
TS4  
[Signature]

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: [;.] GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? No

Test	g/210L	Time
Air Blank	0.000	16:03:45

\*Insufficient Deprivation Period

QAS Card not read, "N" button pressed for Deprivation period TS4

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
11/05/2015

16:04:47

4 of 8  
TS  
[Signature]

---

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:05:39
Diagnostic Test	Pass	16:06:13
Air Blank	0.000	16:06:40
0.100 Cal Check	RFI*	16:06:56
Air Blank	IPS**	16:07:04

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
11/05/2015

16:10:17

50f8  
Tsh  
*[Signature]*

---

Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:11:39
Diagnostic Test	Pass	16:12:13
Air Blank	0.000	16:12:41
0.100 Cal Check	0.000	16:13:01
Air Blank	0.000	16:13:31
Subject Test	0.216*	16:16:39
Air Blank	0.000	16:17:16
Air Blank	0.000	16:17:43
Subject Test	XXX**	16:18:30
Air Blank	0.000	16:19:04

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
11/05/2015 16:19:43

60f8  
B4  
dd

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	16:20:12
0.100 Cal Check	0.098	16:20:32
Air Blank	0.000	16:21:03

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

*Fof8*  
*TD4*  
*Mo*  
*11-25-15*

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 11-5-2015 TIME 1550  
Intoxilyzer Serial # 80-000331 Location DUI-HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C.
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE *Terry Gallegos*

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
11/05/2015 15:50:19

80f8  
T4  
[Signature]

-----  
Standard Lot#: 25414100A4  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:50:49
0.100 Cal Check	0.098	15:51:09
Air Blank	0.000	15:51:39

Preventative Maintenance  
Performed -- Timer Reset