



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

AD  
5-26-16

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5-19-2016 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 4-21-2016 and 5-19-2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on                     .  
The instrument was returned to service on                     .

Terry S. Gallegos

Terry S. Gallegos


5-24-2016

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 24<sup>th</sup>, day of May, 2016.

Elma Escarcega

Notary Public

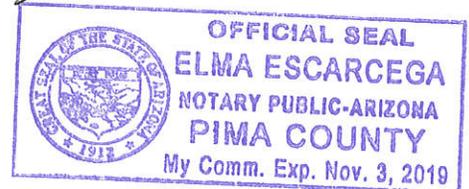

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

10 of 7  
TS4  
5-26-16

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 5-19-2016 TIME 1225  
INTOXILYZER SERIAL # 80-000331 LOCATION DUI-HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
05/19/2016

12:27:04

2057  
704  
JLD

---

Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 12:27:52 |
| Diagnostic Test | Pass   | 12:28:26 |
| Air Blank       | 0.000  | 12:28:53 |
| 0.100 Cal Check | 0.098  | 12:29:14 |
| Air Blank       | 0.000  | 12:29:44 |
| Subject Test    | 0.000  | 12:30:03 |
| Air Blank       | 0.000  | 12:30:31 |
| Wait            | ABT*   |          |
| Air Blank       | 0.000  | 12:31:07 |

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
05/19/2016

12:38:45

30f7  
Tsh  
JL

---

Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 12:40:00 |
| Diagnostic Test | Pass   | 12:40:34 |
| Air Blank       | 0.000  | 12:41:01 |
| 0.100 Cal Check | RFI*   | 12:41:16 |
| Air Blank       | IPS**  | 12:41:26 |

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
05/19/2016

12:42:36

4 of 7  
TJ  
JL

---

Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 12:43:49 |
| Diagnostic Test | Pass   | 12:44:23 |
| Air Blank       | 0.000  | 12:44:51 |
| 0.100 Cal Check | 0.099  | 12:45:12 |
| Air Blank       | 0.000  | 12:45:42 |
| Subject Test    | 0.138* | 12:48:49 |
| Air Blank       | 0.000  | 12:49:22 |
| Air Blank       | 0.000  | 12:49:50 |
| Subject Test    | XXX**  | 12:50:07 |
| Air Blank       | 0.000  | 12:50:42 |

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

5077  
TSG  
JL

INTOXILYZER 8000  
Instrument Initialization  
11:40:47 05/19/2016

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
05/19/2016 12:53:40

-----  
Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 12:54:09 |
| 0.100 Cal Check | 0.098  | 12:54:30 |
| Air Blank       | 0.000  | 12:55:00 |

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6077  
TSB  
ASB  
5-26-16

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 5-19-2016 TIME 1225  
Intoxilyzer Serial # 80000331 Location DUI-HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
05/19/2016 12:25:33

*TG*  
*7 of 7*  
*TG*  
*dd*

-----  
Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 12:26:03 |
| 0.100 Cal Check | 0.098  | 12:26:24 |
| Air Blank       | 0.000  | 12:26:54 |

Preventative Maintenance  
Performed -- Timer Reset