



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

WLL  
8/22/16

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-11-2016 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 7-14-2016 and 8-11-2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Terry Gallegos

Terry S. Gallegos


8-17-2016

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 17<sup>th</sup>, day of August, 2016.

Elma Escarcega

Notary Public

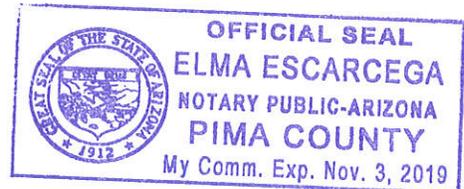

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10f7  
B6  
WLL  
8/22/16

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 8.11.2016 TIME 1020

INTOXILYZER SERIAL # 80-000331 LOCATION DUI-HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE 

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
08/11/2016

10:24:13

2 of 7  
B4  
WLL

---

Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1967  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:25:08
Diagnostic Test	Pass	10:25:43
Air Blank	0.000	10:26:10
0.100 Cal Check	0.098	10:26:31
Air Blank	0.000	10:27:01
Subject Test	0.000	10:27:16
Air Blank	0.000	10:27:43
Wait	ABT*	
Air Blank	0.000	10:28:17

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
08/11/2016

10:29:22

3off  
Tdy  
WLL

Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:30:45
Diagnostic Test	Pass	10:31:19
Air Blank	0.000	10:31:47
0.100 Cal Check	RFI*	10:32:04
Air Blank	IPS**	10:32:18

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
08/11/2016

10:33:13

4 of 7  
TS4  
WLL

Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:34:16
Diagnostic Test	Pass	10:34:50
Air Blank	0.000	10:35:18
0.100 Cal Check	0.098	10:35:38
Air Blank	0.000	10:36:08
Subject Test	0.047*	10:39:16
Air Blank	0.000	10:39:48
Air Blank	0.000	10:40:15
Subject Test	XXX**	10:40:29
Air Blank	0.000	10:41:03

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
08/11/2016 10:41:56

Sof7  
TJ  
WLL

-----  
Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:42:28
0.100 Cal Check	0.097	10:42:48
Air Blank	0.000	10:43:18

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
TS4  
WLL  
8/22/16

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 8-11-2016 TIME 1020  
Intoxilyzer Serial # 80-000331 Location DUI-HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000331  
Core Version: 8105.48  
08/11/2016 10:20:43

*7077  
T54  
WLL*

-----  
Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:21:11
0.100 Cal Check	0.098	10:21:32
Air Blank	0.000	10:22:02

Preventative Maintenance  
Performed -- Timer Reset