



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the
periodic maintenance* for the Intoxilyzer 8000 Serial Number
80-000332 maintained by the Tucson Police Department
Crime Laboratory pursuant to the requirements of the Arizona Department
of Public Safety for a quality assurance program. Tsh

A quality assurance specialist (QAS) tested the instrument on the following
dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2,
for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/22/06 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 12/1/06 and 12/22/06
the instrument was operating properly and accurately.


The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Jennifer L. Coombs
Jennifer L. Coombs

1/4/07

ADPS QAS Permit # 33120

SUBSCRIBED AND SWORN to before me this 4, day of
January, 2007.

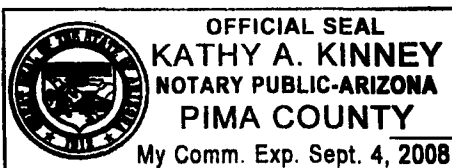

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the
information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Jennifer Coombs 46160 AGENCY Tucson Police Department
DATE December 22, 2006 TIME 1739
INTOXILYZER SERIAL # 80-000332 LOCATION ODE

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed ✓✓
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Jennifer L. Coombs

DPS Form Exh G-4 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

12/22/2006

17:41:36

BH

2 of 8
file

Standard Lot#: 502602I

Last Changed By: J. COOMBS

#46160

QAS: J. COOMBS

#46160

TUCSON PD CRIME LAB

Operator: J. COOMBS

#46160

TUCSON PD CRIME LAB

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:42:42
Diagnostic Test	Pass	17:43:10
Air Blank	0.000	17:43:37
0.100 Cal Check	0.099	17:43:56
Air Blank	0.000	17:44:25
Subject Test	0.000	17:44:41
Air Blank	0.000	17:45:08
Five Minute Wait	ABT*	
Air Blank	0.000	17:45:46

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
12/22/2006

17:46:17

PH

3 of 8
JH

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	17:47:12
Air Blank	0.000	17:47:40

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
12/22/2006

17:48:08

Bh

4 of 8
je

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:49:32
Diagnostic Test	Pass	17:50:00
Air Blank	0.000	17:50:27
0.100 Cal Check	RFI*	17:50:42
Air Blank	0.000	17:51:13

*RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

12/22/2006

17:51:19

BH

5 of 8
ju

Standard Lot#: 502602I

Last Changed By: J. COOMBS

#46160

QAS: J. COOMBS

#46160

TUCSON PD CRIME LAB

Operator: J. COOMBS

#46160

TUCSON PD CRIME LAB

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:52:16
Diagnostic Test	Pass	17:52:43
Air Blank	0.000	17:53:10
0.100 Cal Check	0.099	17:53:29
Air Blank	0.000	17:53:59
Subject Test	0.080*	17:57:07
Air Blank	0.000	17:57:40
Air Blank	0.000	17:58:07
Subject Test	XXX**	17:58:38
Air Blank	0.000	17:59:12

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

12/22/2006

17:59:35

6 of 8
JL
BY


Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	18:00:03
Cal Check	0.098	18:00:22
Air Blank	0.000	18:00:51

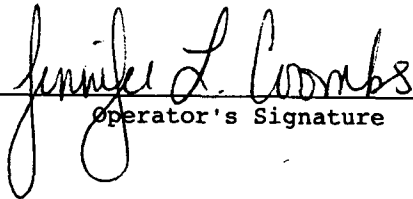

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Jennifer Coombs 46160 AGENCY Tucson Police Department
DATE December 22, 2006 TIME 1739
Intoxilyzer Serial # 80-000332 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Jennifer L. Coombs

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

12/22/2006

17:39:41

Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:40:09
Cal Check	0.099	17:40:27
Air Blank	0.000	17:40:57



Operator's Signature