



CERTIFICATION by CUSTODIAN of RECORDS ^{7/15/07}
Tucson Police Department Crime Laboratory

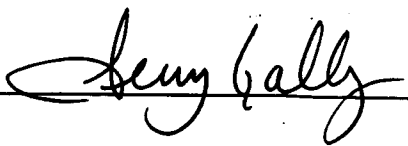
I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5.10.2007 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 4.13.2007 and 5.10.2007 the instrument was operating properly and accurately.

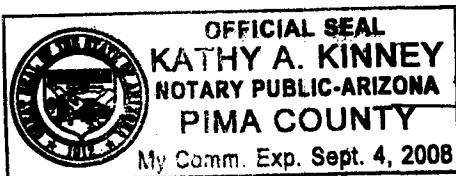
The instrument was taken out of service for repair on .
The instrument was returned to service on .

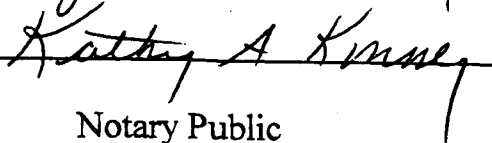

Terry S. Gallegos

5.24.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 24, day of May, 2007.




Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 800010 of 7
TS4^{mm}
5-25-07

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 5.10.2007 TIME 1643
INTOXILYZER SERIAL # 80-000332 LOCATION ODE

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Terry Gallegos

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

05/10/2007

16:44:33

2 of 7
TS4 11/11/11
5-25-07

Standard Lot#: 630601I

Last Changed By: J. COOMBS

#46160

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:45:20
Diagnostic Test	Pass	16:45:47
Air Blank	0.000	16:46:14
0.100 Cal Check	0.101	16:46:33
Air Blank	0.000	16:47:03
Subject Test	0.000	16:47:29
Air Blank	0.000	16:47:57
Five Minute Wait	ABT*	
Air Blank	0.000	16:48:42

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

05/10/2007

16:49:19

3087
TS4⁷⁴₅₋₂₅₆₇

Standard Lot#: 630601I

Last Changed By: J. COOMBS

#46160

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:50:16
Diagnostic Test	Pass	16:50:44
Air Blank	0.000	16:51:11
0.100 Cal Check	RFI*	16:51:29
Air Blank	IPS**	16:51:39

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

05/10/2007

16:51:57

Standard Lot#: 630601I

Last Changed By: J. COOMBS

#46160

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:52:51
Diagnostic Test	Pass	16:53:19
Air Blank	0.000	16:53:46
0.100 Cal Check	0.100	16:54:05
Air Blank	0.000	16:54:34
Subject Test	0.090*	16:57:43
Air Blank	0.000.	16:58:15
Air Blank	0.000	16:58:43
Subject Test	XXX**	16:59:03
Air Blank	0.000	16:59:37

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

40f7
TS4^{MD}₅₋₂₀₀₇

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

05/10/2007

16:59:56

Standard Lot#: 630601I

Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:00:23
Cal Check	0.098	17:00:42
Air Blank	0.000	17:01:11


Operator's Signature

50f7
TSG ¹⁴⁴ 5-25-07

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 80006 of 7
TS4 5/10/07

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 5-10-2007 TIME 1643
Intoxilyzer Serial # 80-000 332 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

05/10/2007

16:43:03

70f7
B45-2501

Standard Lot#: 630601I

Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test g/210L Time

Air Blank 0.000 16:43:31
Cal Check 0.100 16:43:49
Air Blank 0.000 16:44:19


Operator's Signature