



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TJ
6-11-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-8-07 Cal 1 0.099 AC Cal 2 0.100 AC


Between the dates of 5-10-07 and 6-8-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley
Clare L. Conley
6-8-2007

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 8, day of June, 2007.

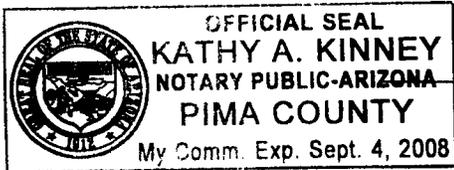

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 8 Cu

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TJ
6-11-07

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley 34517 AGENCY Tucson Police Department
DATE 6-8-07 TIME 0549
INTOXILYZER SERIAL # 80-000332 LOCATION ODG

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Clare Conley

2 of 8 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
06/08/2007 05:51:59

by
6.11.07

Standard Lot#: 630601I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	05:52:52
Diagnostic Test	Pass	05:53:19
Air Blank	0.000	05:53:46
0.100 Cal Check	0.100	05:54:05
Air Blank	0.000	05:54:35
Subject Test	0.000	05:54:53
Air Blank	ABT*	05:55:06
Air Blank	0.000	05:55:34

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 8 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
06/08/2007 05:55:59

Bh
6-11-07

Standard Lot#: 630601I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	05:56:54
Diagnostic Test	Pass	05:57:21
Air Blank	0.000	05:57:48
0.100 Cal Check	0.105	05:58:07
Air Blank	0.000	05:58:37
Subject Test	RFI*	06:01:31
Air Blank	0.000	06:02:06

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 8 cu

TJ
6.11.07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
06/08/2007 06:02:50

Standard Lot#: 630601I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:03:46
Diagnostic Test	Pass	06:04:14
Air Blank	0.000	06:04:41
0.100 Cal Check	0.099	06:05:00
Air Blank	IPS*	06:05:07
Air Blank	0.000	06:05:36

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 8 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
06/08/2007

06:05:57

TJy
6.11.07

Standard Lot#: 630601I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:06:52
Diagnostic Test	Pass	06:07:19
Air Blank	0.000	06:07:46
0.100 Cal Check	0.099	06:08:05
Air Blank	0.000	06:08:35
Subject Test	0.043*	06:11:43
Air Blank	0.000	06:12:16
Air Blank	0.000	06:12:43
Subject Test	XXX**	06:13:00
Air Blank	0.000	06:13:35

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

6 of 8 cu

T84
6-11-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
06/08/2007 06:14:43

Standard Lot#: 630601I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:15:10
Cal Check	0.099	06:15:29
Air Blank	0.000	06:15:59

Clare Conley
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104 (A)

7 of 8 cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TJ
6-11-07

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley 34517 AGENCY Tucson Police Department
DATE 6-8-07 TIME 0549
Intoxilyzer Serial # 80-000332 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

8780a

T54
6-11-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
06/08/2007 05:49:05

Standard Lot#: 630601I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	05:49:32
Cal Check	0.100	05:49:51
Air Blank	0.000	05:50:21

Clare Conley

Operator's Signature