



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TS4
7.10.2007

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80 - 000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-6-07 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 6-8-07 and 7-6-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

7-9-07

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 9, day of July, 2007.

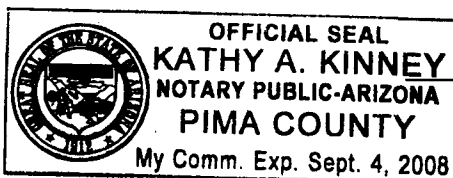

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

B4

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 7-6-07 TIME 0659
INTOXILYZER SERIAL # 80-000332 LOCATION ODE

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS None.

SIGNATURE Clare Conley

2 of 8 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
07/06/2007

07:01:58

TS4

Standard Lot#: 630601I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:02:51
Diagnostic Test	Pass	07:03:18
Air Blank	0.000	07:03:45
0.100 Cal Check	0.099	07:04:04
Air Blank	0.000	07:04:34
Subject Test	0.000	07:04:51
Air Blank	0.000	07:05:18
Five Minute Wait	ABT*	
Air Blank	0.000	07:05:56

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3.18 Cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
07/06/2007

07:06:32

ph

Standard Lot#: 630601I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:07:22
Diagnostic Test	Pass	07:07:50
Air Blank	0.000	07:08:17
0.100 Cal Check	RFI*	07:08:28
Air Blank	0.000	07:08:59

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

07/06/2007

07:09:27

Bq

Standard Lot#: 630601I

Last Changed By: J. COOMBS

#46160

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	07:10:01
Air Blank	0.000	07:10:30

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

548Cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
07/06/2007

07:11:26

ph

Standard Lot#: 630601I
Last Changed By: J. COOMBS #46160
QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SOAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:12:18
Diagnostic Test	Pass	07:12:46
Air Blank	0.000	07:13:13
0.100 Cal Check	0.100	07:13:32
Air Blank	0.000	07:14:02
Subject Test	0.075*	07:17:10
Air Blank	0.000	07:17:43
Air Blank	0.000	07:18:10
Subject Test	XXX**	07:18:29
Air Blank	0.000	07:19:04

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

6 of 8 Cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

07/06/2007

07:19:31

BY

Standard Lot#: 630601I

Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:19:58
Cal Check	0.099	07:20:17
Air Blank	0.000	07:20:47

Clare Conley

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

186

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 7-6-07 TIME 0659
Intoxilyzer Serial # 80-000 332 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$.
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

848Cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

07/06/2007

06:59:20

Standard Lot#: 630601I

Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:59:47
Cal Check	0.100	07:00:06
Air Blank	0.000	07:00:35

Clare Conley

Operator's Signature