



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

Tby
12-21-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/18/07 Cal 10.100 AC Cal 2 0.101 AC

Between the dates of 11/19/07 and 12/18/07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.


Sarah L. Cione

Sarah L. Cione

12/19/07

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 19, day of December, 2007.

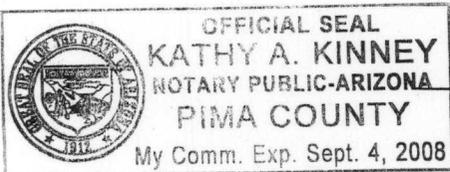

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

1 of 7sc
124

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 12/18/07 TIME 11:25am

INTOXILYZER SERIAL # 80-000332 LOCATION Team 2 - DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
12/18/2007 11:28:42

2 of 7sc
BSH

Standard Lot#: 7254011
Last Changed By: J. OLIVARES #36552

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:29:32
Diagnostic Test	Pass	11:30:01
Air Blank	0.000	11:30:28
0.100 Cal Check	0.101	11:30:47
Air Blank	0.000	11:31:17
Subject Test	0.000	11:31:43
Air Blank	0.000	11:32:10
Five Minute Wait	ABT*	
Air Blank	0.000	11:32:46

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 7 SC
AM

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
12/18/2007 11:35:54

Standard Lot#: 7254011
Last Changed By: J. OLIUARES #36552

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test g/210L Time

Air Blank RFI* 11:36:25
Air Blank IPS** 11:36:33

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 7sc
BY

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
12/18/2007 11:37:19

Standard Lot#: 7254011
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:38:21
Diagnostic Test	Pass	11:38:50
Air Blank	0.000	11:39:17
0.100 Cal Check	0.101	11:39:36
Air Blank	0.000	11:40:06
Subject Test	0.068*	11:43:14
Air Blank	0.000	11:43:47
Air Blank	0.000	11:44:14
Subject Test	XXX**	11:44:34
Air Blank	0.000	11:45:08

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 7 sc

Bh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
12/18/2007 11:45:55

Standard Lot#: 7254011
Last Changed By: J. OLIVARES #36552

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:46:22
Cal Check	0.101	11:46:41
Air Blank	0.000	11:47:11

Sarah Cione
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 75c
B4

QA Specialist Sarah Cione AGENCY Tucson Police Department
 DATE 12/18/07 TIME 11:25am
 Intoxilyzer Serial # 80-000332 Location Team 2 - DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
 OR
 Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 SC
BY

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
12/18/2007 11:25:32

Standard Lot#: 7254011
Last Changed By: J. OLIVARES #36552

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:25:59
Cal Check	0.100	11:26:18
Air Blank	0.000	11:26:48

Sarah Cione

Operator's Signature