



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

Handwritten note: 154 8.31.2009

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 08/20/09 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 07/23/09 and 08/20/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A. The instrument was returned to service on N/A.

Handwritten signature of Cheryl K. Straub

Cheryl K. Straub


Handwritten date: 08/31/09

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 31, day of August, 2009.

Handwritten signature of Notary Public

Notary Public

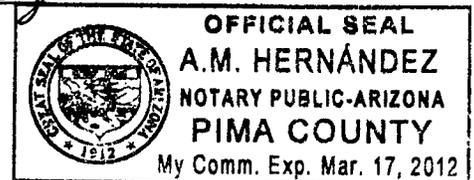

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 8
CB

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

TS4
8-31-09

QA SPECIALIST Cheryl Straub AGENCY Tucson Police Department
DATE 08/20/09 TIME 11:54
INTOXILYZER SERIAL # 80-000332 LOCATION DUT SQUAD HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS RADIO INTERFERENCE OCCURRED DURING DEFICIENT
SAMPLE TEST GIVING A "RFI DETECT." TEST REPEATED.

SIGNATURE Cheryl Straub

288
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
08/20/2009

12:05:44

BY

Standard Lot#: 801002I
Last Changed By: W. HONOMICHL #42992

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:06:50
Diagnostic Test	Pass	12:07:17
Air Blank	0.000	12:07:45
0.100 Cal Check	0.097	12:08:03
Air Blank	0.000	12:08:33
Subject Test	0.000	12:09:22
Air Blank	ABT*	12:09:30
Air Blank	0.000	12:09:58

*Sequence Aborted

Not a Successfully Completed Test Sequence

378
CF

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
08/20/2009

12:11:05

154

Standard Lot#: 801002I
Last Changed By: W. HONOMICHL #42992

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	12:12:02
Air Blank	IPS**	12:12:11

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

498
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
08/20/2009

12:13:10

BY

Standard Lot#: 801002I
Last Changed By: W. HONOMICHL #42992

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:14:10
Diagnostic Test	Pass	12:14:38
Air Blank	0.000	12:15:05
0.100 Cal Check	0.097	12:15:24
Air Blank	0.000	12:15:53
Subject Test	RFI*	12:19:01
Air Blank	0.000	12:19:35

*RFI Detect

Not a Successfully Completed Test Sequence

578
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
08/20/2009

12:20:09

B4

Standard Lot#: 801002I
Last Changed By: W. HONOMICHL #42992

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:21:13
Diagnostic Test	Pass	12:21:41
Air Blank	0.000	12:22:08
0.100 Cal Check	0.097	12:22:27
Air Blank	0.000	12:22:57
Subject Test	0.101*	12:26:05
Air Blank	0.000	12:26:38
Air Blank	0.000	12:27:05
Subject Test	XXX**	12:27:28
Air Blank	0.000	12:28:03

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4
8.31.09

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 08/20/09 TIME 11:54
Intoxilyzer Serial # 80-000332 Location DUT SQUAD HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

172
CG



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

B9
8-31-09

Date: 07-28-09
QAS: William D. Honovich
Intoxilyzer Serial Number: 80-000332
New Lot Number: 801002I
New Expiration Date: 15 Jan 2010

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE [Signature]

282
CS


B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
07/28/2009 18:40:12

Standard Lot#: 8010021
Last Changed By: W. HONOMICHL #42992

OAS: W. HONOMICHL #42992
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	18:40:39
Cal Check	0.098	18:40:58
Air Blank	0.000	18:41:27

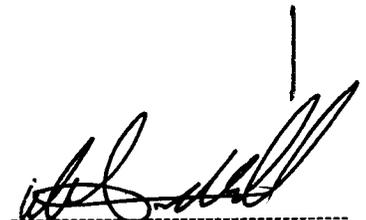

Operator's Signature