



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

SC
6/11/10

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5-26-10 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 4-29-10 and 5-26-10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


6-11-10

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 11, day of June, 2010.

A.M. Hernandez

Notary Public

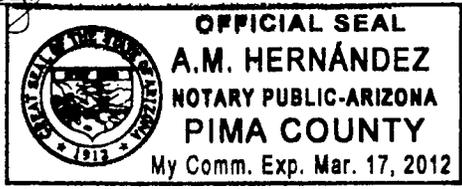

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1587
an
TS4
6.1.2010

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 5-26-10 TIME 0826

INTOXILYZER SERIAL # 80-000332 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
05/26/2010

08:27:56

2007
cu
134

Standard Lot#: 642845
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:28:43
Diagnostic Test	Pass	08:29:10
Air Blank	0.000	08:29:38
0.100 Cal Check	0.098	08:29:56
Air Blank	0.000	08:30:26
Subject Test	0.000	08:30:44
Air Blank	ABT*	08:30:51
Air Blank	0.000	08:31:19

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
05/26/2010 08:31:23

307
cu
154

Standard Lot#: 642845
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:32:12
Diagnostic Test	Pass	08:32:39
Air Blank	0.000	08:33:06
0.100 Cal Check	0.097	08:33:25
Air Blank	RFI*	08:33:35
Air Blank	IPS**	08:33:43

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
05/26/2010

08:34:06

457
cu
BY

Standard Lot#: 642845
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:34:54
Diagnostic Test	Pass	08:35:21
Air Blank	0.000	08:35:48
0.100 Cal Check	0.097	08:36:07
Air Blank	0.000	08:36:37
Subject Test	0.146*	08:39:45
Air Blank	0.000	08:40:18
Air Blank	0.000	08:40:45
Subject Test	XXX**	08:41:17
Air Blank	0.000	08:41:51

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

05/26/2010

08:42:15

Standard Lot#: 642845

Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:42:42
Cal Check	0.097	08:43:01
Air Blank	0.000	08:43:31

Clare Conley

Operator's Signature

557
cu

TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6867
CU
TS4
6.1.2010

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 5-26-10 TIME 0826
Intoxilyzer Serial # 80-000332 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 180 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

05/26/2010

08:26:27

Standard Lot#: 642845

Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:26:54
Cal Check	0.098	08:27:13
Air Blank	0.000	08:27:43

Clare Conley

Operator's Signature

787
cu

T54