



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/22/10 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 6/24/10 and 7/22/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


7/22/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 22, day of July, 20 10.

A.M. Hernandez

Notary Public

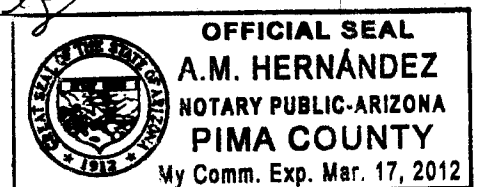

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police DepartmentDATE 7/22/10 TIME 1015INTOXILYZER SERIAL # 80-000332 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE Sarah Cione

DPS Form Exh G-4 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

07/22/2010

10:16:56

2 of 7_{sc}

B4

Standard Lot#: 719919

Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616

TUCSON PD

Operator: S. CIONE #50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:17:41
Diagnostic Test	Pass	10:18:09
Air Blank	0.000	10:18:36
0.100 Cal Check	0.098	10:18:55
Air Blank	0.000	10:19:25
Subject Test	0.000	10:20:03
Air Blank	0.000	10:20:31
Five Minute Wait	ABT*	
Air Blank	0.000	10:21:06

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7_{sc}

134

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

07/22/2010

10:22:28

Standard Lot#: 719919

Last Changed By: W. HONOMICHL

#42992

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:24:06
Diagnostic Test	Pass	10:24:34
Air Blank	0.000	10:25:01
0.100 Cal Check	RFI*	10:25:13
Air Blank	IPS**	10:25:23

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

07/22/2010

10:25:34

4 of 7 SC

B4

Standard Lot#: 719919

Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:26:33
Diagnostic Test	Pass	10:27:00
Air Blank	0.000	10:27:27
0.100 Cal Check	0.097	10:27:46
Air Blank	0.000	10:28:16
Subject Test	0.101*	10:31:24
Air Blank	0.000	10:31:57
Air Blank	0.000	10:32:24
Subject Test	XXX**	10:32:55
Air Blank	0.000	10:33:29

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

07/22/2010

10:34:00

5 of 7 sc

BH

Standard Lot#: 719919

Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:34:27
Cal Check	0.097	10:34:46
Air Blank	0.000	10:35:16

Sarah Cione
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
7-22-2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 7/22/10 TIME 1015
Intoxilyzer Serial # 80-000332 Location DVI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

07/22/2010

10:15:24

7 of 7 SC
BS4

Standard Lot#: 719919

Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:15:51
Cal Check	0.097	10:16:10
Air Blank	0.000	10:16:40

Sarah Cione
Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 2₅

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

TS4
7-22-2010

Date: 7-1-10
QAS: HONOMICHL, WILLIAM
Intoxilyzer Serial Number: 80-000332
New Lot Number: 219919
New Expiration Date: 02-02-2012

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 099 AC.
- (✓) 9. Air blank completed.
- () 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE


2 of 2_{sc}
156

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
07/01/2010 00:05:36

Standard Lot#: 719919
Last Changed By: W. HONOMICHL #42992

OAS: W. HONOMICHL #42992
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	00:06:03
Cal Check	0.099	00:06:22
Air Blank	0.000	00:06:52

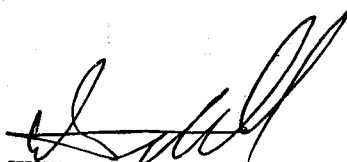
Operator's Signature