



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
4.7.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/31/11 Cal 1 0.098 AC Cal 2 0.096 AC

Between the dates of 3/3/11 and 3/31/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


4/6/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 6th day of April, 20 11.

SMH Hernandez

Notary Public

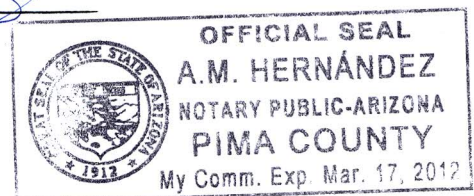

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSY
4-7-2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/31/11 TIME 1156

INTOXILYZER SERIAL # 80-000332 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.096 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

03/31/2011

11:58:01

BM

Standard Lot#: 719919

Last Changed By: W. HONOMICHL

#42992

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:58:43
Diagnostic Test	Pass	11:59:11
Air Blank	0.000	11:59:38
0.100 Cal Check	0.097	11:59:57
Air Blank	0.000	12:00:26
Subject Test	0.000	12:00:49
Air Blank	0.000	12:01:16
Five Minute Wait	ABT*	
Air Blank	0.000	12:02:06

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

03/31/2011

12:03:28

754

Standard Lot#: 719919

Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616

TUCSON PD

Operator: S. CIONE #50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:04:15
Diagnostic Test	Pass	12:04:43
Air Blank	RFI*	12:04:53
Air Blank	IPS**	12:05:01

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

03/31/2011

12:05:26

Standard Lot#: 719919

Last Changed By: W. HONOMICHL

#42992

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:06:08
Diagnostic Test	Pass	12:06:35
Air Blank	0.000	12:07:02
0.100 Cal Check	0.097	12:07:21
Air Blank	0.000	12:07:51
Subject Test	0.048*	12:10:59
Air Blank	0.000	12:11:31
Air Blank	0.000	12:11:58
Subject Test	XXX**	12:12:18
Air Blank	0.000	12:12:51

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

03/31/2011

12:13:35

Standard Lot#: 719919

Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:14:02
Cal Check	0.096	12:14:21
Air Blank	0.000	12:14:51

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

TS4
4.7.2011

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/31/11 TIME 1156
Intoxilyzer Serial # 80-000332 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

03/31/2011

11:56:31

126

Standard Lot#: 719919

Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:56:58
Cal Check	0.098	11:57:17
Air Blank	0.000	11:57:47

Sarah Cione

Operator's Signature