



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
11-7-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/17/11 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 9/15/11 and 10/17/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


10/27/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27, day of October, 2011.

Ant Hernandez

Notary Public

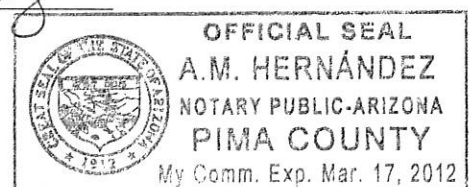

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

BS
11.7.2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 10/17/11 TIME 1041
INTOXILYZER SERIAL # 80-000332 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS 1 day past due on Calibration Check

SIGNATURE Sarah Cione

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

10/17/2011

10:46:58

Standard Lot#: 04911100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:47:56
Diagnostic Test	Pass	10:48:24
Air Blank	0.000	10:48:51
0.100 Cal Check	0.097	10:49:10
Air Blank	0.000	10:49:39
Subject Test	0.000	10:49:54
Air Blank	0.000	10:50:21
Five Minute Wait	ABT*	
Air Blank	0.000	10:51:15

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 9_{sc}

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

10/17/2011

10:51:32

B4

Standard Lot#: 04911100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:52:06
Air Blank	IPS**	10:52:14

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

4 of 9_{sc}

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

10/17/2011

10:53:02

TS4

Standard Lot#: 04911100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:53:43
Diagnostic Test	Pass	10:54:11
Air Blank	0.000	10:54:38
0.100 Cal Check	0.097	10:54:57
Air Blank	0.000	10:55:26
Subject Test	0.074*	10:58:35
Air Blank	0.000	10:59:07
Air Blank	0.000	10:59:34
Subject Test	XXX**	11:00:18
Air Blank	0.000	11:00:52

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 9_{sc}

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

10/17/2011

11:01:16

TSy

Standard Lot#: 04911100A1

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:01:44
Cal Check	0.097	11:02:02
Air Blank	0.000	11:02:32

Sarah Cione

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/17/11 TIME 1041
Intoxilyzer Serial # 80-000332 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 9_{sc}

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

10/17/2011

10:41:54

TS4

Standard Lot#: 719919

Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:42:22
Cal Check	0.097	10:42:41
Air Blank	0.000	10:43:10

Sarah Cione

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

8 of 9sc

TS4
11.7-2011

Date: 10/17/11

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-000332

New Lot Number: 04911100A1

New Expiration Date: 3/1/2013

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0.098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
10/17/2011

10:45:18

TS4

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:45:45
Cal Check	0.098	10:46:04
Air Blank	0.000	10:46:34

Sarah Cione

Operator's Signature