



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
 7.15.2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/27/13 Cal 1 0.097 AC Cal 2 0.095 AC

Between the dates of 5/30/13 and 6/27/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

7/11/13


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 11, day of July, 2013.

A.M. Hernandez

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

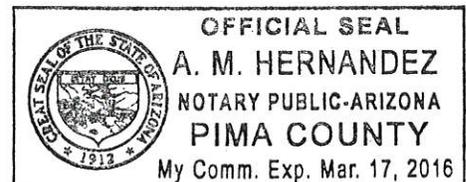

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7 sc  
TS4  
7-15-2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/27/13 TIME 903

INTOXILYZER SERIAL # 80-000332 LOCATION ODS - Tm 1

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 095 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
06/27/2013

09:05:31

2 of 7<sub>sc</sub>  
TS4

---

Standard Lot#: 13712100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:06:14
Diagnostic Test	Pass	09:06:49
Air Blank	0.000	09:07:16
0.100 Cal Check	0.097	09:07:37
Air Blank	0.000	09:08:07
Subject Test	0.000	09:08:29
Air Blank	0.000	09:08:57
Wait	ABT*	
Air Blank	0.000	09:09:34

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
06/27/2013

09:09:37

3 of 7 SC  
Tsh

---

Standard Lot#: 13712100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:10:21
Diagnostic Test	Pass	09:10:56
Air Blank	0.000	09:11:24
0.100 Cal Check	0.096	09:11:45
Air Blank	0.000	09:12:14
Subject Test	RFI*	09:12:27
Air Blank	IPS**	09:12:35

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
06/27/2013

09:12:39

4 of 7  
TS4

---

Standard Lot#: 13712100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:13:20
Diagnostic Test	Pass	09:13:55
Air Blank	0.000	09:14:23
0.100 Cal Check	0.097	09:14:43
Air Blank	0.000	09:15:13
Subject Test	0.445*	09:18:21
Air Blank	0.000	09:18:54
Air Blank	0.000	09:19:21
Subject Test	XXX**	09:19:37
Air Blank	0.000	09:20:12

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
06/27/2013 09:20:22

5 of 7  
SC  
TBH

-----  
Standard Lot#: 13712100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:20:50
0.100 Cal Check	0.095	09:21:11
Air Blank	0.000	09:21:41

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

6 of 7 SC  
TS4  
7.15.2013

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 6/27/13 TIME 903  
Intoxilyzer Serial # 80-000332 Location ODS-Tm 1

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
06/27/2013 09:03:49

7 of 7  
5c  
TSh

-----  
Standard Lot#: 13712100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:04:18
0.100 Cal Check	0.097	09:04:39
Air Blank	0.000	09:05:08

Preventative Maintenance  
Performed -- Timer Reset