



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
8/8/13

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7.25.2013 Cal 1 0.095 AC Cal 2 0.094 AC

Between the dates of 6-27-2013 and 7-25-2013 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

Terry Gallegos

Terry S. Gallegos


8-7-2013

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 7<sup>th</sup>, day of August, 2013.

Elma Escarcega

Notary Public

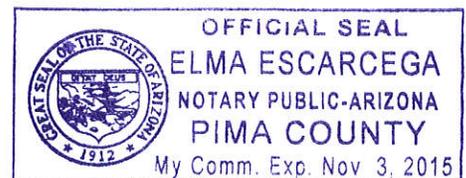

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 9  
TJG

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
8/8/13

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 7-25-2013 TIME 1240

INTOXILYZER SERIAL # 80-000332 LOCATION ODS

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 094 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/25/2013

12:42:23

2 of 9  
T04  
SC

---

Standard Lot#: 13712100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:43:12
Diagnostic Test	Pass	12:43:47
Air Blank	0.000	12:44:15
0.100 Cal Check	0.095	12:44:35
Air Blank	0.000	12:45:05
Subject Test	0.000	12:45:22
Air Blank	0.000	12:45:49
Wait	ABT*	
Air Blank	0.000	12:46:28

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/25/2013

12:46:39

30f9  
TSh  
SC

---

Standard Lot#: 13712100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP

DOB: 11/11/1987

Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:47:27
Diagnostic Test	Pass	12:48:05
Air Blank	0.000	12:48:32
0.100 Cal Check	RFI*	12:48:48
Air Blank	IPS**	12:48:58

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/25/2013

12:50:23

4 of 9  
TJH  
SC

---

Standard Lot#: 13712100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:51:14
Diagnostic Test	Pass	12:51:49
Air Blank	0.000	12:52:17
0.100 Cal Check	0.096	12:52:38
Air Blank	0.000	12:53:07
Subject Test	0.375*	12:56:15
Air Blank	0.000	12:56:49
Air Blank	0.000	12:57:16
Subject Test	XXX**	12:57:33
Air Blank	0.000	12:58:08

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/25/2013 12:58:21

50f9  
TD4

-----  
Standard Lot#: 13712100A2  
Last Changed By: T. GALLEGOS #28847  
-----

SC

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:58:51
0.100 Cal Check	0.094	12:59:12
Air Blank	0.000	12:59:41

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

60f 9  
TJL  
sc  
8/8/13

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 7-25-2013 TIME 1240  
Intoxilyzer Serial # 80-000332 Location ODS

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 095 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/25/2013 12:40:42

70f9  
TS4

-----  
Standard Lot#: 13712100A2  
Last Changed By: T. GALLEGOS #28847  
-----

5c

-----  
QAS: S. CIONE #50616  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:41:12
0.100 Cal Check	0.095	12:41:33
Air Blank	0.000	12:42:03

Preventative Maintenance NOT  
Performed -- Timer NOT Reset



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

80f9

TJH

SC  
8/8/13

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

Date: 7.25.2013

QAS: Terry Gallegos 28847

Intoxilyzer Serial Number: 80-000332

New Lot Number: 304 12100A1

New Expiration Date: 12.1.2014

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
  - Enter Target Value: 0.100
  - Enter Control Lot #: From the tank label
  - Enter Exp Date: From the tank label (mm/dd/yyyy)
  - Edit info Y/N? Y or N
  - Update Std Info? Y ENTER
  - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 096 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/25/2013 13:05:38

9 of 9  
TS4  
SC  
8/8/13

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:06:09
0.100 Cal Check	0.096	13:06:30
Air Blank	0.000	13:07:00

Preventative Maintenance  
Performed -- Timer Reset