



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

5C
9/16/13

I hereby certify that these documents are a true and correct copy of the
periodic maintenance\* for the Intoxilyzer 8000 Serial Number
80-000332 maintained by the Tucson Police Department
Crime Laboratory pursuant to the requirements of the Arizona Department
of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following
dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2,
for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-22-2013 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 7-25-2013 and 8-22-2013
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on
The instrument was returned to service on

[Handwritten signature of Terry S. Gallegos]

Terry S. Gallegos


9.10.2013

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 10, day of
September, 2013.

[Handwritten signature of A.M. Hernandez]

Notary Public

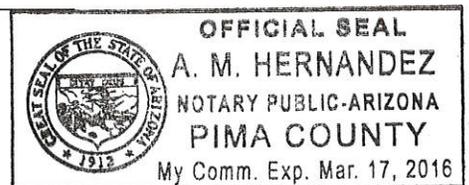

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1057  
TSh

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
9/16/13

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 8-22-2013 TIME 1216

INTOXILYZER SERIAL # 80-000.332 LOCATION ODS

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 0.096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS \_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
08/22/2013

12:18:05

20f7  
DH SC

---

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:19:00
Diagnostic Test	Pass	12:19:35
Air Blank	0.000	12:20:02
0.100 Cal Check	0.098	12:20:23
Air Blank	0.000	12:20:53
Subject Test	0.000	12:21:24
Air Blank	0.000	12:21:51
Wait	ABT*	
Air Blank	0.000	12:22:38

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
08/22/2013

12:22:42

30f7  
TS4 SC

---

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:23:45
Diagnostic Test	Pass	12:24:20
Air Blank	0.000	12:24:48
0.100 Cal Check	RFI*	12:25:02
Air Blank	IPS**	12:25:12

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
08/22/2013

12:25:17

4 of 7  
TDH  
SC

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TPD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:26:20
Diagnostic Test	Pass	12:26:55
Air Blank	0.000	12:27:22
0.100 Cal Check	0.097	12:27:43
Air Blank	0.000	12:28:13
Subject Test	0.088*	12:31:20
Air Blank	0.000	12:31:54
Air Blank	0.000	12:32:21
Subject Test	XXX**	12:32:42
Air Blank	0.000	12:33:17

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
08/22/2013 12:33:30

*Sof 7*  
*TS4*  
*SC*

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:34:00
0.100 Cal Check	0.096	12:34:21
Air Blank	0.000	12:34:51

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

60 of 7  
TSh

SC  
9/16/13

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 8.22.2013 TIME 1216  
Intoxilyzer Serial # 80-000332 Location ODS

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
08/22/2013 12:16:20

7057  
Dh

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

50

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	12:16:48
0.100 Cal Check	0.097	12:17:08
Air Blank	0.000	12:17:38