



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

SL
10/2/13

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9-18-2013 Cal 1 0.097 AC Cal 2 0.094 AC

Between the dates of 8-22-2013 and 9-18-2013 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.


Terry S. Gallegos

9-26-2013

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 26th, day of September, 2013.

Notary Public

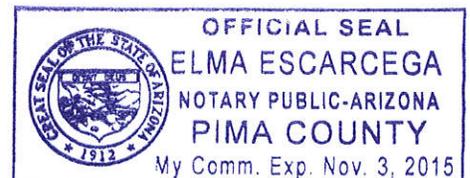

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10 of 7
By SC
10/2/13

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 9.18.2013 TIME 1205

INTOXILYZER SERIAL # 80-000332 LOCATION ODS

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 094 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
09/18/2013

12:08:16

2077
PL SC

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 12:09:17 |
| Diagnostic Test | Pass | 12:09:52 |
| Air Blank | 0.000 | 12:10:19 |
| 0.100 Cal Check | 0.097 | 12:10:40 |
| Air Blank | 0.000 | 12:11:10 |
| Subject Test | 0.000 | 12:11:31 |
| Air Blank | 0.000 | 12:11:58 |
| Wait | ABT* | |
| Air Blank | 0.000 | 12:12:36 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
09/18/2013

12:12:41

3087
TS4

5C

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 12:13:37 |
| Diagnostic Test | Pass | 12:14:12 |
| Air Blank | 0.000 | 12:14:39 |
| 0.100 Cal Check | RFI* | 12:14:56 |
| Air Blank | IPS** | 12:15:04 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
09/18/2013

12:15:17

4 of 7
TS4
SC

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 12:16:07 |
| Diagnostic Test | Pass | 12:16:42 |
| Air Blank | 0.000 | 12:17:09 |
| 0.100 Cal Check | 0.097 | 12:17:30 |
| Air Blank | 0.000 | 12:18:00 |
| Subject Test | 0.259* | 12:21:07 |
| Air Blank | 0.000 | 12:21:41 |
| Air Blank | 0.000 | 12:22:08 |
| Subject Test | XXX** | 12:22:32 |
| Air Blank | 0.000 | 12:23:07 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
09/18/2013 12:23:22

507
T84
5c

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 12:23:52 |
| 0.100 Cal Check | 0.094 | 12:24:13 |
| Air Blank | 0.000 | 12:24:42 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6047
TS4
SC
10/2/13

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 9.18.2013 TIME 1205
Intoxilyzer Serial # 80-000332 Location ODS

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
09/18/2013 12:05:16

7077
TS4
SC

Standard Lot#: 30412100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 12:05:43 |
| 0.100 Cal Check | 0.097 | 12:06:04 |
| Air Blank | 0.000 | 12:06:34 |