



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10-17-2013 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 9-18-2013 and 10-17-2013 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

[Handwritten signature of Terry S. Gallegos]

Terry S. Gallegos


11-7-2013

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 7, day of November, 2013.

[Handwritten signature of A. M. Hernandez]

Notary Public

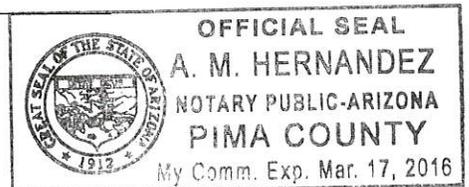

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
TS4 SC  
11/19/13

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 10-17-2013 TIME 1155  
INTOXILYZER SERIAL # 80-000332 LOCATION ODS

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
10/17/2013

11:57:21

20f7  
TS4 SC

---

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:58:13
Diagnostic Test	Pass	11:58:48
Air Blank	0.000	11:59:15
0.100 Cal Check	0.098	11:59:36
Air Blank	0.000	12:00:06
Subject Test	0.000	12:00:37
Air Blank	0.000	12:01:04
Wait	ABT*	
Air Blank	0.000	12:01:46

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
10/17/2013

12:01:56

3077  
TS4 SC

---

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:02:46
Diagnostic Test	Pass	12:03:21
Air Blank	0.000	12:03:48
0.100 Cal Check	RFI*	12:04:04
Air Blank	IPS**	12:04:14

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
10/17/2013

12:04:25

4 of 7  
TS4 SC

---

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:05:13
Diagnostic Test	Pass	12:05:48
Air Blank	0.000	12:06:15
0.100 Cal Check	0.097	12:06:36
Air Blank	0.000	12:07:06
Subject Test	0.120*	12:10:13
Air Blank	0.000	12:10:48
Air Blank	0.000	12:11:15
Subject Test	XXX**	12:11:40
Air Blank	0.000	12:12:14

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
10/17/2013 12:15:55

507  
TJH SC

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:16:26
0.100 Cal Check	0.096	12:16:46
Air Blank	0.000	12:17:16

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6077  
TS4  
SC  
11/19/13

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 10-17-2013 TIME 1155  
Intoxilyzer Serial # 80-000332 Location ODS

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
10/17/2013

11:55:19

7 of 7  
Dh

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

SC

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	11:55:47
0.100 Cal Check	0.097	11:56:07
Air Blank	0.000	11:56:37