



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

5C  
8/12/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7.24.2014 Cal 1 0.095 AC Cal 2 0.093 AC

Between the dates of 6.25.2014 and 7.24.2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 7.24.2014.  
The instrument was returned to service on 7.24.2014.  
Cal 3 0.100 AC

Terry Gallegos

Terry S. Gallegos


8.12.2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 12, day of August, 20 14.

A.M. Hernandez

Notary Public

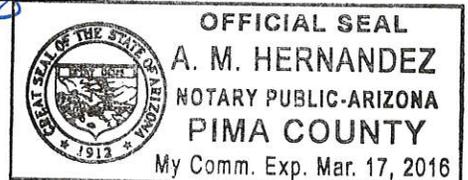

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10F12  
TS4  
SC  
8/12/14

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 7-24-2014 TIME 1329  
INTOXILYZER SERIAL # 80-000332 LOCATION ODS

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 093 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/24/2014

13:30:45

20812  
TS4  
SC

---

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:31:35
Diagnostic Test	Pass	13:32:10
Air Blank	0.000	13:32:38
0.100 Cal Check	0.094	13:32:59
Air Blank	0.000	13:33:28
Subject Test	0.000	13:34:00
Air Blank	0.000	13:34:28
Wait	ABT*	
Air Blank	0.000	13:35:02

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/24/2014

13:35:06

3 of 12  
TSh

SC

---

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:35:54
Diagnostic Test	Pass	13:36:29
Air Blank	0.000	13:36:56
0.100 Cal Check	RFI*	13:37:11
Air Blank	IPS**	13:37:20

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/24/2014

13:37:22

4 of 12  
TS4

---

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

SC

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:38:10
Diagnostic Test	Pass	13:38:45
Air Blank	0.000	13:39:13
0.100 Cal Check	0.094	13:39:33
Air Blank	0.000	13:40:03
Subject Test	0.454*	13:43:11
Air Blank	0.000	13:43:45
Air Blank	0.000	13:44:12
Subject Test	XXX**	13:45:00
Air Blank	0.000	13:45:34

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/24/2014 13:45:46

*Sof 12*  
*TS4*  
*SC*

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:46:17
0.100 Cal Check	0.093	13:46:38
Air Blank	0.000	13:47:08

Preventative Maintenance  
Performed -- Timer Reset

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/24/2014 13:51:01

60712  
TS4

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

SC

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:51:32
0.100 Cal Check	0.093	13:51:53
Air Blank	0.000	13:52:22

Preventative Maintenance  
Performed -- Timer Reset

Removed and reinstalled gas tank TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

70512  
TS4

3C  
8/12/14

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 7.24.2014 TIME 1329  
Intoxilyzer Serial # 80-000332 Location ODS

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0.095 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/24/2014 13:29:11

8 of 12  
TS4  
SC

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	13:29:39
0.100 Cal Check	0.095	13:30:00
Air Blank	0.000	13:30:29



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

9 of 12  
TS4

INTOXILYZER SERVICE RECORD

SC  
8/2/14

Date: 7.24.2014

QAS: Terry Gallegos 28847

Intoxilyzer Serial Number: 80-000332

Reason For Service: Calibration checks are low  
Brought to Crime Lab for service

Date Out of Service: 7.24.2014

Calibration Check completed prior to removal from service

Date: 7.24.2014

SQAP completed prior to removal from service

Date: 7.24.2014

Corrective Action:

Sent to CMI, Inc.

Replaced regulator O-ring and gas tank change.  
Returned to service at ODS.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/24/2014 17:10:09

10 of 12  
TS4

SC

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	17:10:39
0.100 Cal Check	0.094	17:11:00
Air Blank	0.000	17:11:30

Preventative Maintenance NOT  
Performed -- Timer NOT Reset

Changed O-ring on regulator  
TS4



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

11 of 12  
TS4

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

SC  
8/12/14

Date: 7.24.2014  
QAS: Terry Gallegos 28847  
Intoxilyzer Serial Number: 80-000332  
New Lot Number: 11514100A2  
New Expiration Date: 05.01.2016

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
  - Enter Target Value: 0.100
  - Enter Control Lot #: From the tank label
  - Enter Exp Date: From the tank label (mm/dd/yyyy)
  - Edit info Y/N? Y or N
  - Update Std Info? Y ENTER
  - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/24/2014 17:14:15

12 of 12  
TS4

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

SC

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	17:14:45
0.100 Cal Check	0.100	17:15:05
Air Blank	0.000	17:15:35

Preventative Maintenance  
Performed -- Timer Reset

Changed tank TS4