



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

154
2.10.2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/5/2015 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 1/7/2015 and 2/5/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis


2/5/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 5th, day of February, 2015.

Elma Escarcega

Notary Public

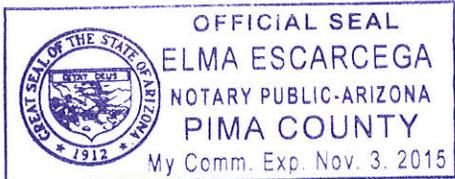

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
WLL
TBH
2.10.2015

QA SPECIALIST WILL LOOMIS 53545 AGENCY Tucson Police Department

DATE 2/5/2015 TIME 10:04

INTOXILYZER SERIAL # 80-000332 LOCATION ODS

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
02/05/2015

10:07:14

2 of 7
WLL
TS4

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:08:12 |
| Diagnostic Test | Pass | 10:08:47 |
| Air Blank | 0.000 | 10:09:14 |
| 0.100 Cal Check | 0.100 | 10:09:35 |
| Air Blank | 0.000 | 10:10:05 |
| Subject Test | 0.000 | 10:10:20 |
| Air Blank | 0.000 | 10:10:48 |
| Wait | ABT* | |
| Air Blank | 0.000 | 10:11:23 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
02/05/2015

10:12:03

3 of 7
WLL

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:13:01 |
| Diagnostic Test | Pass | 10:13:36 |
| Air Blank | 0.000 | 10:14:04 |
| 0.100 Cal Check | RFI* | 10:14:18 |
| Air Blank | IPS** | 10:14:26 |

*RFI Detect

**Improper Sample

TS4

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
02/05/2015

10:14:45

4 of 7
WLL
TS4

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:15:38 |
| Diagnostic Test | Pass | 10:16:13 |
| Air Blank | 0.000 | 10:16:41 |
| 0.100 Cal Check | 0.100 | 10:17:02 |
| Air Blank | 0.000 | 10:17:31 |
| Subject Test | 0.593* | 10:20:39 |
| Air Blank | 0.000 | 10:21:13 |
| Air Blank | 0.000 | 10:21:40 |
| Subject Test | XXX** | 10:21:54 |
| Air Blank | 0.000 | 10:22:28 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
02/05/2015 10:23:20

5 of 7
WLV
134

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:23:51 |
| 0.100 Cal Check | 0.098 | 10:24:12 |
| Air Blank | 0.000 | 10:24:42 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
WLL
TS4
2.10.2015

QA Specialist WILL LOOMIS 53345 AGENCY Tucson Police Department
DATE 2/5/2015 TIME 10:04
Intoxilyzer Serial # 80-000332 Location ODS

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
02/05/2015

10:04:51

7 of 7
wll
154

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:05:27 |
| 0.100 Cal Check | 0.100 | 10:05:48 |
| Air Blank | 0.000 | 10:06:18 |

Preventative Maintenance
Performed -- Timer Reset