



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

Handwritten initials and date:
 10/5/15

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9.18.2015 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 8.21.2015 and 9.18.2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
 The instrument was returned to service on .

Handwritten signature of Terry S. Gallegos

Terry S. Gallegos


9.29.2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 29th, day of September, 2015.

Handwritten signature of Elma Escarcega

Notary Public

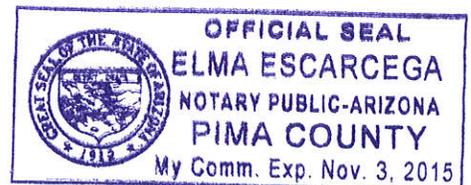

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8
TS4
DLS
10/5/15

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 9.18.2015 TIME 1724

INTOXILYZER SERIAL # 80-000332 LOCATION ODS

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Terry Gallegos

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
09/18/2015

17:25:44

20f8
154
dlb

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 17:27:02 |

*Sequence Aborted

30f8
TS4
de

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
09/18/2015

17:27:07

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 17:27:53 |
| Diagnostic Test | Pass | 17:28:28 |
| Air Blank | 0.000 | 17:28:55 |
| 0.100 Cal Check | 0.099 | 17:29:16 |
| Air Blank | 0.000 | 17:29:46 |
| Subject Test | 0.000 | 17:30:00 |
| Air Blank | 0.000 | 17:30:28 |
| Wait | ABT* | |
| Air Blank | 0.000 | 17:31:02 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
09/18/2015

17:31:05

4 of 8
TS4
JW

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 17:31:56 |
| Diagnostic Test | Pass | 17:32:31 |
| Air Blank | 0.000 | 17:32:58 |
| 0.100 Cal Check | RFI* | 17:33:13 |
| Air Blank | IPS** | 17:33:20 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Sof8
TS4
JL

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
09/18/2015 17:36:51

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 17:37:36 |
| Diagnostic Test | Pass | 17:38:11 |
| Air Blank | 0.000 | 17:38:38 |
| 0.100 Cal Check | 0.098 | 17:38:59 |
| Air Blank | 0.000 | 17:39:28 |
| Subject Test | 0.191* | 17:42:36 |
| Air Blank | 0.000 | 17:43:09 |
| Air Blank | 0.000 | 17:43:36 |
| Subject Test | XXX** | 17:44:16 |
| Air Blank | 0.000 | 17:44:51 |

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Core Version: 8105.48
Serial Number: 80-000332
09/18/2015

17:24:00

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 17:24:44 |
| 0.100 Cal Check | 0.099 | 17:25:04 |
| Air Blank | 0.000 | 17:25:34 |
| End of test | | |

6 of 8
TS4
JLW

Printer error, reprinted calibration report
TS4

Preventive Maintenance NOT
Performed -- Timer NOT Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8
TS4
Add
10/5/15

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 9-18-2015 TIME 1724
Intoxilyzer Serial # 80-000332 Location ODS

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
09/18/2015 17:24:00

80f8
B4
[Signature]

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 17:24:44 |
| 0.100 Cal Check | 0.099 | 17:25:04 |
| Air Blank | 0.000 | 17:25:34 |

Preventative Maintenance
Performed -- Timer Reset