



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the <sup>PM</sup> periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.


Date: 12/1/06 Cal 1 0.097 AC Cal 2 0.099 AC

The instrument was placed into service on 12/1/06.

Jennifer L. Coombs  
Jennifer L. Coombs  
12/27/06

ADPS QAS Permit # 33120

SUBSCRIBED AND SWORN to before me this 27, day of December, 2006.

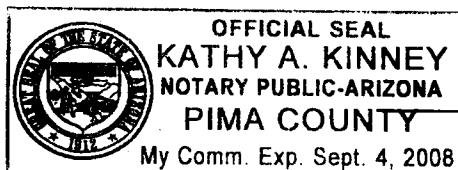

Kathy A. Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000STANDARD QUALITY ASSURANCE PROCEDURE *for*QA SPECIALIST Jennifer Coombs 461160 AGENCY Tucson Police DepartmentDATE December 1, 2006 TIME 1225INTOXILYZER SERIAL # 80000334 LOCATION LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- 
- (✓) 2. Date check.

## OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
- 
- (✓) 2. Error recognition logic system functioning.
- 
- Not a Successfully Completed Test Sequence printed ✓✓
- 
- (✓) 3. Proper sample recognition system.
- 
- Not a Successfully Completed Test Sequence printed ✓
- 
- Deficient sample printed. ✓
- 
- (✓) 4. Standard Calibration Check standard 0.100 AC
- 
- Result: 0.097 AC

Instrument is operating properly and accurately. YES ✓ NO       

## COMMENTS

---

---

---

---

SIGNATURE

Jennifer L. Coombs

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
12/01/2006

12:34:00

164

2 of 8  
JH

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160  
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160  
TUCSON PD CRIME LAB

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:36:00
Diagnostic Test	Pass	12:36:27
Air Blank	0.000	12:36:55
0.100 Cal Check	0.098	12:37:13
Air Blank	0.000	12:37:44
Subject Test	0.000	12:38:18
Air Blank	0.000	12:38:45
Five Minute Wait	ABT*	
Air Blank	0.000	12:39:54

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000334

12/01/2006

12:40:00

Bh

3 of 8  
file

Standard Lot#: 502602I

Last Changed By: J. COOMBS

#46160

QAS: J. COOMBS

#46160

TUCSON PD CRIME LAB

Operator: J. COOMBS

#46160

TUCSON PD CRIME LAB

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	12:41:13
Air Blank	0.000	12:41:42

\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000334

12/01/2006

12:42:08

B4

4 of 8  
file

Standard Lot#: 502602I

Last Changed By: J. COOMBS

#46160

QAS: J. COOMBS

#46160

TUCSON PD CRIME LAB

Operator: J. COOMBS

#46160

TUCSON PD CRIME LAB

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:43:02
Diagnostic Test	Pass	12:43:29
Air Blank	0.000	12:43:57
0.100 Cal Check	RFI*	12:44:11
Air Blank	0.000	12:44:43

\*RFI Detect

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000334

12/01/2006

12:45:21

BH

5 of 8  
file

Standard Lot#: 502602I

Last Changed By: J. COOMBS

#46160

QAS: J. COOMBS

#46160

TUCSON PD CRIME LAB

Operator: J. COOMBS

#46160

TUCSON PD CRIME LAB

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:46:25
Diagnostic Test	Pass	12:46:53
Air Blank	0.000	12:47:20
0.100 Cal Check	0.099	12:47:39
Air Blank	0.000	12:48:09
Subject Test	0.045*	12:51:18
Air Blank	0.000	12:51:50
Air Blank	0.000	12:52:17
Subject Test	XXX**	12:52:57
Air Blank	0.000	12:53:31

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

12/01/2006

12:55:43

6 of 8  
JL

-----  
Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

PM

-----  
QAS: J. COOMBS #46160

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	12:56:10
Cal Check	0.097	12:56:29
Air Blank	0.000	12:56:59

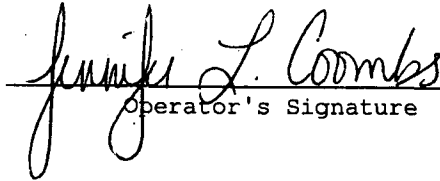
  
\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

7 of 8  
JL

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES *for*  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Jennifer Coombs 46160 AGENCY Tucson Police Department  
DATE December 1, 2006 TIME 1225  
Intoxilyzer Serial # 80000334 Location LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Jennifer L. Coombs

DPS Form Exh G-2 (Rev. 05-01)



INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

12/01/2006

12:25:46


Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	12:26:13
Cal Check	0.099	12:26:32
Air Blank	0.000	12:27:02

  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

TM

1 of 1  
file

Intoxilyzer Serial Number: 80-000334

Initial Service Date: 12-1-06

CMI, Inc. Certificate of Calibration Date: 11-6-03, 5-2-05