



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

*[Signature]*

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 90-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-28-2006 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 12-1-2006 and 12-28-2006 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.


*[Signature: Terry S. Gallegos]*

Terry S. Gallegos

1-4-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 4, day of January, 2007.

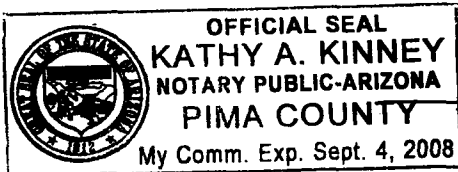

*[Signature: Kathy A. Kinney]*  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 8  
P4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 2002 AGENCY Tucson Police DepartmentDATE 12-28-2006 TIME 1704INTOXILYZER SERIAL # 80-000334 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

## OPERATIONAL TESTS


- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed ✓  
Deficient sample printed. ✓  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ✓ NO    

## COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

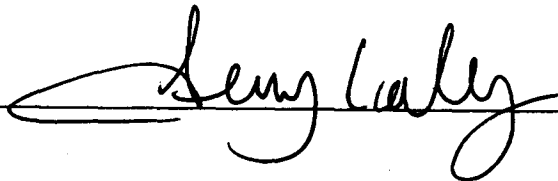

2 of 8  
Pg

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
12/28/2006

17:06:20

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:07:59
Diagnostic Test	Pass	17:08:27
Air Blank	0.000	17:08:54
0.100 Cal Check	0.098	17:09:13
Air Blank	0.000	17:09:43
Subject Test	0.000	17:09:59
Air Blank	0.000	17:10:26
Five Minute Wait	ABT*	
Air Blank	0.000	17:11:09

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

30f8  
TJ4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
12/28/2006 17:11:20

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:12:14
Diagnostic Test	Pass	17:12:41
Air Blank	0.000	17:13:08
0.100 Cal Check	0.098	17:13:27
Air Blank	0.000	17:13:57
Reference	IPS*	17:14:04
Air Blank	0.000	17:14:32

\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
12/28/2006

17:15:58

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	17:16:36
Air Blank	0.000	17:17:04

\*RFI Detect

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
12/28/2006

17:17:24

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:19:01
Diagnostic Test	Pass	17:19:29
Air Blank	0.000	17:19:56
0.100 Cal Check	0.098	17:20:15
Air Blank	0.000	17:20:45
Subject Test	0.094*	17:23:54
Air Blank	0.000	17:24:26
Air Blank	0.000	17:24:53
Subject Test	XXX**	17:25:31
Air Blank	0.000	17:26:06

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

12/28/2006

17:26:45

60f8  
B4

Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:27:13
Cal Check	0.097	17:27:31
Air Blank	0.000	17:28:02

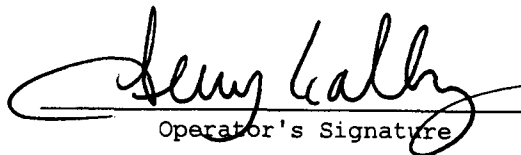
  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104 (A)

70F8  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

jk

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 20847 AGENCY Tucson Police Department  
DATE 12-28-2006 TIME 1704  
Intoxilyzer Serial # 80-000334 Location DVI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)



INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

12/28/2006

17:04:38

Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:05:05
Cal Check	0.098	17:05:24
Air Blank	0.000	17:05:54

  
Operator's Signature