



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

7/14/07  
650

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000 334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5.24.2007 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 4.26.2007 and 5.24.2007 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

Terry S. Gallegos

6.5.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 5, day of June, 2007.

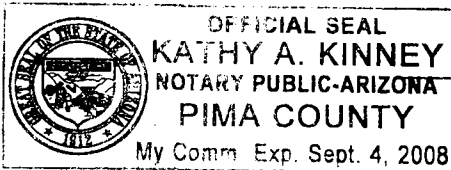

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 80001 of 7  
Dy 6-5-07

## STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 5.24.2007 TIME 1630  
INTOXILYZER SERIAL # 80-000334 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

## OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO    

## COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Terry Gallegos

DPS Form Exh G-4 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000334

07/24/1987

16:33:52

20fz  
Tdy 6:567

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS

#28847

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Name	g.c.c.	Time
Air Blank	0.000	16:34:41
Diagnostic Test	Pass	16:35:08
Air Blank	0.000	16:35:35
0.100 Cal Check	0.099	16:35:54
Air Blank	0.000	16:36:24
Subject Test	0.000	16:36:54
15 Minute Wait	ABP	16:37:22
Air Blank	0.000	16:38:18

^Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

1111 TUCSON POLICE DEPT

Officer Number: 80000334

09/24/2007

16:38:30

3087  
TSG  
6-5-07

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS

#28847

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1977

Sex: F

Weight: 111

Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:39:23
Diagnostic Test	Pass	16:39:51
Air Blank	RFI*	16:40:05
Sample	TP**	16:40:15

RFI Detect:

\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INNOVIZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000334

05/24/2007

16:40:22

4 of 7  
Tdy  
6-5-07

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS

#28847

Q101 T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:41:16
Air Blank	0.000	16:41:43
Air Blank	0.000	16:42:11
0.100 Cal. Check	0.099	16:42:29
Air Blank	0.000	16:43:00
Subject Test	0.075*	16:46:08
Air Blank	0.000	16:46:40
Air Blank	0.000	16:47:07
Subject Test	XXX**	16:47:28
Air Blank	0.000	16:48:01

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

05/24/2007

16:48:31

Soft  
By 6-5-07

Operator: T. GALLEGOS #28847

Q28: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/mL	Time
1st Blank	0.000	16:48:33
2nd Blank	0.000	16:48:33
Air Blank	0.000	16:49:47

  
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

## ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

## STANDARD CALIBRATION CHECK PROCEDURE

6077  
TSY 6507QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 5.24.2007 TIME 1630  
Intoxilyzer Serial # 80-000334 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8100

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

05/04/2007

16:30:58

70F7  
TS4  
2/23  
6-5-07

Standard Lot#: 6297000

Test Charged Port: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/dL	Time
Air Blank	0.000	16:31:25
		16:31:42
		16:32:14

  
Operator's Signature