



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
1-2-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/30/08 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 12/18/08 and 12/30/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 12/30/08.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


1/2/2009

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 2, day of January, 2009.

A.M. Hernandez

Notary Public

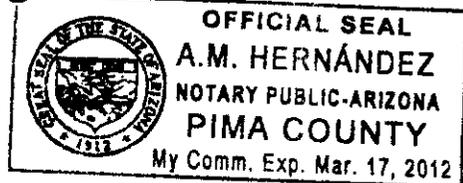

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 9 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

T54
1-2-09

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 12/30/08 TIME 13:32
INTOXILYZER SERIAL # 80-000334 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000334
12/30/2008

13:34:02

2 of 9 SC
B4

Standard Lot#: 811201I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:34:55
Diagnostic Test	Pass	13:35:22
Air Blank	0.000	13:35:50
0.100 Cal Check	0.100	13:36:08
Air Blank	0.000	13:36:38
Subject Test	0.000	13:37:09
Air Blank	0.000	13:37:36
Five Minute Wait	ABT*	
Air Blank	0.000	13:38:23

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000334
12/30/2008

13:38:32

3 of 9 sc
p4

Standard Lot#: 811201I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	13:39:07
Air Blank	IPS**	13:39:15

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000334
12/30/2008

13:39:30

4 of 9 SC
B4

Standard Lot#: 811201I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:40:15
Diagnostic Test	Pass	13:40:43
Air Blank	0.000	13:41:10
0.100 Cal Check	0.100	13:41:29
Air Blank	0.000	13:41:58
Subject Test	0.042*	13:45:07
Air Blank	0.000	13:45:38
Air Blank	0.000	13:46:06
Subject Test	XXX**	13:46:38
Air Blank	0.000	13:47:11

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
12/30/2008 13:48:09

5 of 9 SC
B4

Standard Lot#: 811201I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:48:37
Cal Check	0.100	13:48:56
Air Blank	0.000	13:49:25

Smah Cione
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 9 sc

TS4
1-2-09

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 12/30/08 TIME 13:32
Intoxilyzer Serial # 80-000334 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
12/30/2008

13:32:12

7 of 9 SC

184

Standard Lot#: 811201I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:32:40
Cal Check	0.101	13:32:58
Air Blank	0.000	13:33:28

Sarah Cione

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9 sc

INTOXILYZER SERVICE RECORD

TS4
1-2-09

Date: 12/30/08

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-000334

Reason For Service: Officer notified us that instrument was failing the modern check. At the lab, we could not perform a subject test as when "y" was entered for evidentiary test, instrument would display "Tank below minimum pressure" when tank was full & connected. Replaced regulator for SQAP. Modern test passed diagnostics at lab, just had problem w/ regulator.

Date Out of Service: 12/30/08

Calibration Check completed prior to removal from service

Date: 12/30/08

SQAP completed prior to removal from service

Date: 12/30/08

Corrective Action:

Sent to CMI, Inc.

Data was downloaded



9 of 9 sc

104

Return Material Authorization Form

This form MUST be completed and enclosed with item(s) being returned.
Failure to complete and return this form may cause delays in crediting / repairing.

1. Item(s) Returned: Instrument Supplies Other
(Note: please ship items in their original shipping container(s) or a similar protective box.)

2. Instrument Model: Intox 8000 Serial Number 80-000334

3. Bill To Address: Finance / Accounts Payable
P.O. Box 27450
TUCSON, AZ 85726-7450
Ship To Address: (NO PO BOXES)
Tucson Police Department
Crime Lab
270 S. Stone Ave
TUCSON, AZ 85701-1917

Customer # 105075
Name Terry Gallegos
Phone: (520)-791-4494 Fax: (520) 791-4642
Email: Terry.Gallegos@tucsonaz.gov

4. Reason for Return: Note - If you are sending an item for repair, please give a detailed description of the problem. Please list any special instructions that you may have concerning this return. Note - Returns for credit may be subject to a 15% restocking charge.
Originally failed moderm test multiple times. Also regulator / regulator connector not recognized, as could not perform an evidentiary test. Instrument would flag low tank pressure. Another regulator worked fine.

5. Purchase Order Number (attach a copy of P.O. if applicable) 20669

6. I Authorize All Repairs: Yes No

7. I Authorize Repairs Up To: \$ _____

8. I Require An Estimate Regardless of Cost Yes No

Note: An evaluation fee will apply to estimates that are not repaired.

9. Authorized By: Terry Gallegos
Name (Please Print)

[Signature] Date 12.31.08
Signature