



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

BSG  
6-26-09

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/24/09 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of N/A and N/A the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on 6/24/09.

Sarah Cione

Sarah L. Cione


6/24/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 24<sup>th</sup>, day of June, 20 09.

A.M. Hernandez

Notary Public

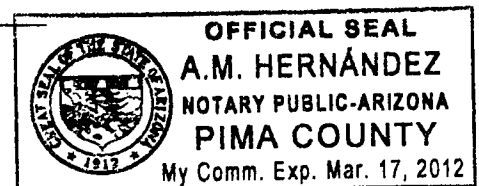

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 12 sc

## ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

BSH  
6.26.09

## STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department  
DATE 6/24/09 TIME 14:59  
INTOXILYZER SERIAL # 80-000334 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

## OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO       

## COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Sarah Cione

2 of 12<sub>sc</sub>

184

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000334

06/24/2009

15:01:59

---

Standard Lot#: 826001I

Last Changed By: C. CONLEY

#34517

---

QAS: S. CIONE

#50616

TUCSON PD

---

Operator: S. CIONE

#50616

TUCSON PD

---

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	15:02:56
Diagnostic Test	Pass	15:03:24
Air Blank	0.000	15:03:51
0.100 Cal Check	0.099	15:04:10
Air Blank	0.000	15:04:40
Subject Test	0.000	15:05:14
Air Blank	0.000	15:05:41
Five Minute Wait	ABT*	
Air Blank	0.000	15:06:17

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 12 sc

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
06/24/2009

15:06:41

134

Standard Lot#: 826001I  
Last Changed By: C. CQNLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	15:07:17
Air Blank	IPS**	15:07:26

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 12 sc  
B4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
06/24/2009

15:07:41

---

Standard Lot#: 826001I  
Last Changed By: C. CONLEY #34517

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	15:08:26
Diagnostic Test	Pass	15:08:54
Air Blank	0.000	15:09:21
0.100 Cal Check	0.099	15:09:40
Air Blank	0.000	15:10:10
Subject Test	0.075*	15:13:18
Air Blank	0.000	15:13:50
Air Blank	0.000	15:14:17
Subject Test	XXX**	15:14:48
Air Blank	0.000	15:15:22

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

06/24/2009

15:15:50

-----  
Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616.

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:16:18
Cal Check	0.099	15:16:36
Air Blank	0.000	15:17:06

Sarah Cione

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 12<sub>sc</sub>

## ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000154  
6-26-09

## STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 6/24/09 TIME 14:59  
Intoxilyzer Serial # 80-000334 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

06/24/2009

14:59:16


-----  
Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616.

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:59:43
Cal Check	0.100	15:00:02
Air Blank	0.000	15:00:32

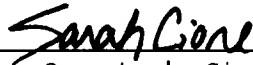

Operator's Signature





TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

8 of 12<sub>sc</sub>

INTOXILYZER SERVICE RECORD

154  
6-26-09

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-000334

Date Returned to Service: 6/24/09



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

9 of 12  
5c

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

Date: 6-18-09  
QAS: Clare Conley  
Intoxilyzer Serial Number: 80-000334  
New Lot Number: 826001I  
New Expiration Date: 9/24/10

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:
  - Enter Target Value: 0.100
  - Enter Control Lot #: From the tank label
  - Enter Exp Date: From the tank label (mm/dd/yyyy)
  - Edit info Y/N? Y or N
  - Update Std Info? Y ENTER
  - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 099 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

10 of 12sc

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

06/18/2009

10:39:59

-----  
Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:40:26
Cal Check	0.099	10:40:45
Air Blank	0.000	10:41:15

  
\_\_\_\_\_  
Operator's Signature

319-7026 BH

B  
I  
L  
L  
  
T  
OSHIP  
T O

MODEL # 8000	S/N 80 000334	DATE RCVD.	P.O. # 10667	REQ.	BILL CODE 1	SERVICE LOCATION CODE
PROD. CODE	EXTRA PARTS RCVD. F. 10112					

DESCRIPTION OF PROBLEM MOTOR / TRANSMISSION REVERSE FAILURES	<input type="checkbox"/> QUOTE	WARRANTY <input type="checkbox"/> SERVICE DEPT. <input type="checkbox"/> NEW
---	--------------------------------	---

MODEM / RELAY / PULSE SOURCE / PULSE RATE

☒ PRELIMINARY TEST      ☒ DISMANTLE, INSPECT, DIAGNOSTIC TEST  
POWER SUPPLY      ☒ TEST, INSPECT      ☐ REPLACE

☐ MAIN PWB's - REMOVE & CLEAN, REINSTALL, TEST  
PWB's      ☒ DIAGNOSE, REPAIR      ☐ REPLACE

☒ PWB's - SURFACE CLEAN IN PLACE

☐ FILTER WHEEL - REMOVE, CLEAN, REINSTALL, CHECK CELL INTERIOR  
CHOPPER WHEEL MOTOR      ☐ TEST & INSPECT      ☐ REPLACE

☒ CELL - REMOVE, CLEAN, REINSTALL, REPLACE LENSES/GASKETS REQ'D  
IR LIGHT SOURCE      ☐ REMOVE, REALIGN      ☐ REPLACE  
PRINTER ASSY.      ☐ DISASSEMBLE      ☐ CLEAN & LUBRICATE  
                                 ☐ REMOVE & REPLACE BELT      ☐ REALIGNMENT  
                                 ☐ REPLACE      ☒ TEST PRINTER FUNCTIONS

BATTERIES - REPLACE      ☐ CPU      ☐ RAM

☐ REPLACE PARTS      ☒ REASSEMBLE      ☒ CALIBRATION      ☒ FINAL TEST

☐ MISC. \_\_\_\_\_

PBT      ☐ PRELIMINARY TESTING - PRESSURE, BATTERY, ZERO SET  
                 ☐ DISASSEMBLE      ☐ REPLACE PARTS      ☐ REASSEMBLE  
                 ☐ TEST      ☐ CALIBRATE      ☐ FINAL TEST

TOXITEST      ☐ INCOMING LEAK TEST      ☐ VISUAL/MECH. INSPECTION  
                 ☐ TEST THERMOMETER & THERMOSTAT      ☐ DISASSEMBLE  
                 ☐ REPLACE PARTS      ☐ REASSEMBLE      ☐ OUTGOING & LEAK TEST  
                 ☐ FINAL TEST

☐ MISC. \_\_\_\_\_

[illegible]**COMMENTS:**

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

**HOURS:**

TECH.

DATE \_\_\_\_\_

**CAL TECH.**

DATE \_\_\_\_\_

FORM # CMI 8WO 1 REV (-)  
WHITE-SERVICE DEPT/YELLOW-PACKING SLIP/PINK-INVOICE COPY

**SHIPPING DATE/CHARGES**

SPECIAL SHIPPING INSTRUCTIONS

☐ COD      ☐ AMOUNT \$

# CERTIFICATE OF CALIBRATION

## INTOXILYZER® 8000

*This is to certify that the calibration of INTOXILYZER serial number 80-000334, manufactured by CMI, subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with the National Safety Administration Standard for Devices to Measure Breath Alcohol ( F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R126 Edition 1998 (E). Calibration solutions are traceable to NIST (NBS) standard material 1828.*

Date 06/3/09

Signed Kristi Robinson  
Technician



316 EAST 9TH STREET  
OWENSBORO, KY. 42303