



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
5.31.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/26/11 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 4/28/11 and 5/26/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -

The instrument was returned to service on -

Sarah Cione


Sarah L. Cione

5/30/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 30, day of

MAY, 20 11

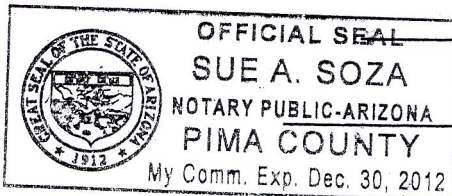

Sue A Soza

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

TS4
5.31.2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 5/26/11 TIME 1232
INTOXILYZER SERIAL # 80-000334 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7_{5c}

154

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000334

05/26/2011

12:34:03

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:34:48
Diagnostic Test	Pass	12:35:16
Air Blank	0.000	12:35:43
0.100 Cal Check	0.100	12:36:02
Air Blank	0.000	12:36:31
Subject Test	0.000	12:37:06
Air Blank	0.000	12:37:34
Five Minute Wait	ABT*	
Air Blank	0.000	12:38:46

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7sc

134

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000334
05/26/2011 12:40:38

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	12:41:19
Air Blank	IPS**	12:41:26

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7_{SC}

154

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000334

05/26/2011

12:41:33

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:42:16
Diagnostic Test	Pass	12:42:44
Air Blank	0.000	12:43:11
0.100 Cal Check	0.099	12:43:30
Air Blank	0.000	12:44:00
Subject Test	0.050*	12:47:08
Air Blank	0.000	12:47:39
Air Blank	0.000	12:48:07
Subject Test	XXX**	12:48:45
Air Blank	0.000	12:49:18

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

05/26/2011

12:50:53

Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:51:20
Cal Check	0.099	12:51:39
Air Blank	0.000	12:52:09

Sarah Cione

Operator's Signature

5 of 7_{sc}

134

EXHIBIT G-2

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

TS4
5-31-2011

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 5/26/11 TIME 1232
Intoxilyzer Serial # 80-000334 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution-0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

05/26/2011

12:32:24

Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:32:51
Cal Check	0.099	12:33:10
Air Blank	0.000	12:33:39

Sarah Cione

Operator's Signature