



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

5C
2/14/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-6-2014 Cal 1 0.099 AC Cal 2 0.097 AC

Between the dates of 1-9-2014 and 2-6-2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry Gallegos

Terry S. Gallegos

2-13-2014


ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 13th, day of February, 20 14.

Elma Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

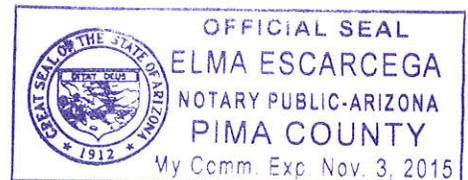

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10 of 7
TS4
sc
2/14/14

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 2.6.2014 TIME 1121

INTOXILYZER SERIAL # 80-000334 LOCATION Sgts Office-ODM

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7
TS4
SC

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
02/06/2014 11:22:37

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

OAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SOAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:23:28
Diagnostic Test	Pass	11:24:04
Air Blank	0.000	11:24:32
0.100 Cal Check	0.099	11:24:52
Air Blank	0.000	11:25:22
Subject Test	0.000	11:25:52
Air Blank	0.000	11:26:20
Wait	ABT*	
Air Blank	0.000	11:27:04

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3087
TSh
SC

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
02/06/2014 11:33:37

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

DAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SOAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:34:27
Diagnostic Test	Pass	11:35:03
Air Blank	0.000	11:35:30
0.100 Cal Check	0.099	11:35:51
Air Blank	RFI*	11:35:58
Air Blank	IPS**	11:36:06

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4087
754

SC

INTOXILYZER #000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
02/06/2014 11:36:20

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TPD

Subject: TEST, SOAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:37:34
Diagnostic Test	Pass	11:38:11
Air Blank	0.000	11:38:38
0.100 Cal Check	0.099	11:38:59
Air Blank	0.000	11:39:28
Subject Test	0.178*	11:42:36
Air Blank	0.000	11:43:09
Air Blank	0.000	11:43:36
Subject Test	XXX**	11:44:23
Air Blank	0.000	11:44:57

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully
Completed
Test Sequence

50f 7
TS4
sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
02/06/2014 11:45:34

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:46:04
0.100 Cal Check	0.097	11:46:25
Air Blank	0.000	11:46:54

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

60 of 7
TS4
SC
2/14/14

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 2-6-2014 TIME 1121
Intoxilyzer Serial # 80-000334 Location Sgts Office - ODM

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

7 of 7
TD
SC

INTOXILYZER 8000
Instrument Initialization
10:47:27 02/06/2014

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
02/06/2014 11:21:02

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

DAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	11:21:30
0.100 Cal Check	0.099	11:21:50
Air Blank	0.000	11:22:20