



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TSh  
 4.10.2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/2/2015 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 3/5/2015 and 4/2/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on                     .  
 The instrument was returned to service on                     .

William L Loomis

William L. Loomis


4/2/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 2nd, day of April, 2015.

Elma Escarcega

Notary Public

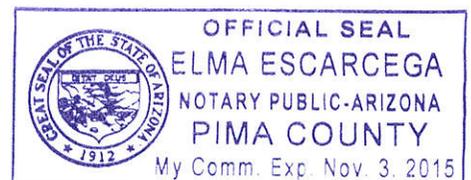

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7 will  
Bk 4.10.2015

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 4/2/2015 TIME 08:40

INTOXILYZER SERIAL # 80-000334 LOCATION ODM - SGT'S OFFICE

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Will Loomis

INTOX ILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000334  
Core Version: 8105.48  
04/02/2015 08:43:16

2 of 7  
WLL  
TS4

Standard Lot#: 11514100A2  
Last Changed By: T. CALLEGOS #28847

OAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SOAP, TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:44:37
Diagnostic Test	Pass	08:45:13
Air Blank	0.000	08:45:41
0.100 Cal Check	0.098	08:46:02
Air Blank	0.000	08:46:31
Subject Test	0.000	08:46:47
Air Blank	0.000	08:47:14
Wait	ABT*	
Air Blank	0.000	08:48:08

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000334  
Core Version: 8105.48  
04/02/2015 08:50:50

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:52:14
Diagnostic Test	Pass	08:52:51
Air Blank	0.000	08:53:18
0.100 Cal Check	RFI*	08:53:31
Air Blank	IPS**	08:53:45

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

3 of 7  
WLL  
TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000334  
Core Version: 8105.48  
04/02/2015 08:54:52

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

OAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SOAP, TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:55:49
Diagnostic Test	Pass	08:56:25
Air Blank	0.000	08:56:53
0.100 Cal Check	0.097	08:57:13
Air Blank	0.000	08:57:43
Subject Test	0.000*	09:00:50
Air Blank	0.000	09:01:24
Air Blank	0.000	09:01:52
Subject Test	XXX**	09:02:16
Air Blank	0.000	09:02:50

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

4 of 7  
wll  
TS4

5 of 7  
WLL  
TSG

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000334  
Core Version: 8105.48  
04/02/2015 09:03:52

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

OAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:04:22
0.100 Cal Check	0.096	09:04:43
Air Blank	0.000	09:05:13

Preventative Maintenance  
Performed -- Timer Reset

]

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

6 of 7 will  
TS4 4.10.2015

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 4/2/2015 TIME 08:40  
Intoxilyzer Serial # 80-000334 Location ODM-5GT'S OFFICE

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8011  
Location: TUCSON POLICE DEPT  
Serial Number: 81-000334  
Core Version: 8105.48  
04/02/2015 08:40:29

7 of 7 will  
TSH

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:41:00
0.100 Cal Check	0.097	08:41:20
Air Blank	0.000	08:41:50

Preventative Maintenance  
Performed -- Timer Reset