



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

*Wm*

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1-26-2007 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 12-28-2006 and 1-26-2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.


*Terry S. Gallegos*

Terry S. Gallegos

2-27-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 27, day of February, 2007.

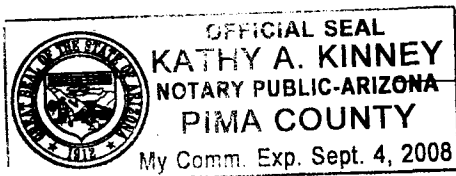

*Kathy A. Kinney*  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

10f 9  
B4 m

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 1-26-2007 TIME 1601  
INTOXILYZER SERIAL # 80-000712 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

( ✓ ) 1. Clock time check.  
( ✓ ) 2. Date check.

OPERATIONAL TESTS

( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.  
( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓  
( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ✓ NO       

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
01/26/2007

16:03:15

20f 9  
154

---

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:04:08
Diagnostic Test	Pass	16:04:36
Air Blank	0.000	16:05:03
0.100 Cal Check	0.097	16:05:22
Air Blank	0.000	16:05:51
Subject Test	0.000	16:06:24
Air Blank	0.000	16:06:51
Five Minute Wait	ABT*	
Air Blank	IPS**	16:08:17

\*Sequence Aborted

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
01/26/2007

16:09:24

3 of 9  
B4

---

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111

Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:10:11
Diagnostic Test	Pass	16:10:39
Air Blank	0.000	16:11:06
0.100 Cal Check	RFI*	16:11:21
Air Blank	0.000	16:11:52

\*RFI Detect

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
01/26/2007

16:12:28

4059  
TS4

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:13:19
Diagnostic Test	Pass	16:13:46
Air Blank	0.000	16:14:14
0.100 Cal Check	0.098	16:14:33
Air Blank	0.000	16:15:02
Subject Test	0.044*	16:18:10
Air Blank	0.000	16:18:41
Air Blank	0.000	16:19:08
Subject Test	XXX**	16:19:25
Air Blank	0.000	16:19:59

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000712

01/26/2007

16:20:44

-----  
Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160  
-----

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	16:21:11
Cal Check	0.097	16:21:30
Air Blank	0.000	16:22:00

  
-----  
Operator's Signature

Sof 9  
TSC

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

6 of 9  
pg 1

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 1-26-2007 TIME 1601  
Intoxilyzer Serial # 80-000712 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000712

01/26/2007

16:01:37

7 of 9  
TS6

-----  
Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160  
-----

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	16:02:04
Cal Check	0.098	16:02:23
Air Blank	0.000	16:02:53

  
-----  
Operator's Signature





TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

8 of 9  
Dy

Date: 1.26.2007  
QAS: Terry Gallegos 28847  
Intoxilyzer Serial Number: 80-000712  
New Lot Number: 5075016  
New Expiration Date: 3.18.2007

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000712

01/26/2007

16:29:34

9 of 9  
TS4

-----  
Standard Lot#: 507501G

Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	16:30:02
Cal Check	0.098	16:30:20
Air Blank	0.000	16:30:50

  
Operator's Signature