



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

MM  
6-5-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5.24.2007 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 4.26.2007 and 5.24.2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on         .  
The instrument was returned to service on         .


Terry S. Gallegos

Terry S. Gallegos

6.5.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 5, day of June, 2007.

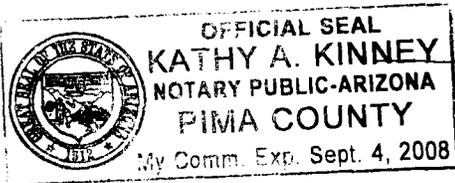

Kathy A Kinney

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 800010 of 7  
TS4  
6/20/07

## STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
 DATE 5-24-2007 TIME 1602  
 INTOXILYZER SERIAL # 80-000712 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
 (✓) 2. Date check.

## OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
 (✓) 2. Error recognition logic system functioning.  
 Not a Successfully Completed Test Sequence printed ✓✓  
 (✓) 3. Proper sample recognition system.  
 Not a Successfully Completed Test Sequence printed  
 Deficient sample printed.  
 (✓) 4. Standard Calibration Check standard 0. 100 AC  
 Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO 

## COMMENTS

Officer Galvan reported that the space bar on the  
keyboard is warped. Will order keyboard replacement  
from CMI, Inc.

## SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
05/24/2007

16:04:13

2087  
TSY  
6-5-07

---

Standard Lot#: 635301I  
Last Changed By: G. EPPLEY #26227

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:05:15
Diagnostic Test	Pass	16:05:42
Air Blank	0.000	16:06:09
0.100 Cal Check	0.101	16:06:28
Air Blank	0.000	16:06:58
Subject Test	0.000	16:07:25
Air Blank	0.000	16:07:52
Five Minute Wait	ABT*	
Air Blank	0.000	16:08:27

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
05/24/2007

16:08:38

3 of 7  
74  
6/5/07

---

Standard Lot#: 635301I  
Last Changed By: G. EPPLEY #26227

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:09:28
Diagnostic Test	Pass	16:09:55
Air Blank	0.000	16:10:22
0.100 Cal Check	0.100	16:10:41
Air Blank	0.000	16:11:10
Subject Test	RFI*	16:11:21
Air Blank	IPS**	16:11:31

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
05/24/2007 16:11:42

4 of 7  
TSY  
6-307

Standard Lot#: 635301I  
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP,  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:12:34
Diagnostic Test	Pass	16:13:01
Air Blank	0.000	16:13:29
0.100 Cal Check	0.101	16:13:47
Air Blank	0.000	16:14:17
Subject Test	0.036*	16:17:25
Air Blank	0.000	16:17:56
Air Blank	0.000	16:18:24
Subject Test	XXX**	16:19:14
Air Blank	0.000	16:19:48

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.



THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

*6 of 7  
TSJ  
5-5-07*

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 5.24.2007 TIME 1602  
Intoxilyzer Serial # 80-000712 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE *Terry Gallegos*

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000712  
05/24/2007 16:02:22

7 of 7  
TS4  
6-5-07

-----  
Standard Lot#: 635301I  
Last Changed By: G. EPPLEY #26227  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	16:02:50
Cal Check	0.100	16:03:08
Air Blank	0.000	16:03:38

  
-----  
Operator's Signature