



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TS4
6-12-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-7-07 Cal 1 0.101 AC Cal 2 0.102 AC

Between the dates of 5-24-07 and 6-7-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

6-8-2007

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 8, day of June, 2007.

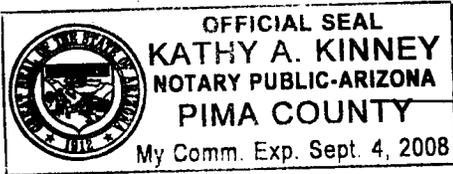

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

~~1 of 8 cu cu~~

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104 (A)

1 of 9 cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

TSK
6.12.07

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley 34517 AGENCY Tucson Police Department

DATE 6-7-07 TIME 1152

INTOXILYZER SERIAL # 80-000712 LOCATION Traffic

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Clare Conley

289 u
B4
6-12-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
06/07/2007 11:05:50

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:06:52
Diagnostic Test	Pass	11:07:19
Air Blank	0.000	11:07:47
0.100 Cal Check	0.100	11:08:05
Air Blank	0.000	11:08:35
Subject Test	XXX*	11:08:52
Air Blank	0.000	11:09:21

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 9 cu

B4
6-12-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
06/07/2007 11:09:44

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:10:43
Diagnostic Test	Pass	11:11:10
Air Blank	0.000	11:11:38
0.100 Cal Check	0.100	11:11:57
Air Blank	0.000	11:12:26
Subject Test	0.000	11:12:47
Air Blank	ABT*	11:12:58
Air Blank	0.000	11:13:27

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 9 Cu

TJG
6-12-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
06/07/2007 11:13:41

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:14:49
Diagnostic Test	Pass	11:15:16
Air Blank	0.000	11:15:44
0.100 Cal Check	0.100	11:16:02
Air Blank	IPS*	11:16:11
Air Blank	0.000	11:16:39

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 9 cu

TS4
6-12-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
06/07/2007 11:17:16

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:18:20
Diagnostic Test	Pass	11:18:47
Air Blank	0.000	11:19:14
0.100 Cal Check	RFI*	11:19:26
Air Blank	0.000	11:19:56

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

659 cu

T34
6-12-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
06/07/2007 11:20:09

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:21:09
Diagnostic Test	Pass	11:21:36
Air Blank	0.000	11:22:03
0.100 Cal Check	0.100	11:22:22
Air Blank	0.000	11:22:52
Subject Test	0.020	11:23:07
Air Blank	0.000	11:23:38
Five Minute Wait		
Air Blank	0.000	11:28:35
Subject Test	0.000	11:28:57
Air Blank	0.000	11:29:24
Five Minute Wait		
Air Blank	0.000	11:34:24
Subject Test	0.073*	11:37:33
Air Blank	0.000	11:38:05

*Deficient Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7069 cu

TS4
6-12-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
06/07/2007 11:38:55

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:39:23
Cal Check	0.100	11:39:41
Air Blank	0.000	11:40:11

Clare Conley
Operator's Signature

8 of 9 cu

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

B4
6-12-07

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley 34517 AGENCY Tucson Police Department
DATE 6-7-07 TIME 1152
Intoxilyzer Serial # 80-000712 Location Traffic

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

9 of 9 cu

By
6-12-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
06/07/2007 11:04:10

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:04:37
Cal Check	0.100	11:04:56
Air Blank	0.000	11:05:25

Clare Conley

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

CU
6-12-07

10/1
754

Date: 5.30.2007

QAS: Terry Gallegos 28847

Intoxilyzer Serial Number: 80-000712

Reason For Service: The space bar key on the keyboard is misshapen from the heat of the light used by the officer. I replaced the keyboard with a spare keyboard.

Date Out of Service: N/A

Calibration Check completed prior to removal from service

Date: _____

SQAP completed prior to removal from service

Date: _____

Corrective Action: Sent to Applied Electronics

