



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

pu
10-25-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9-5-2007 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 8-2-2007 and 9-5-2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .


Terry S. Gallegos

Terry S. Gallegos

10-25-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 25, day of October, 2007.

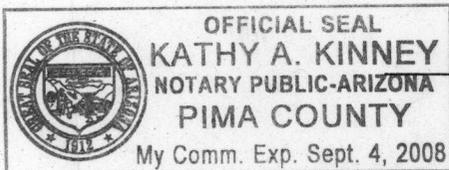

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

CU
10.25.07
10 of 7
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 9-5-2007 TIME 1722

INTOXILYZER SERIAL # 80-000712 LOCATION Lab

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 090 AC

Instrument is operating properly and accurately. YES NO

COMMENTS The instrument was picked up from DUI HQ
on 9-5-2007 for SQAP.
Ex G-2 was done at 34 days.

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
09/05/2007

17:25:40

CW
2 of 7
BY

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:26:32
Diagnostic Test	Pass	17:26:59
Air Blank	0.000	17:27:27
0.100 Cal Check	0.098	17:27:45
Air Blank	0.000	17:28:15
Subject Test	0.000	17:28:41
Air Blank	0.000	17:29:09
Five Minute Wait	ABT*	
Air Blank	0.000	17:29:44

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 16-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
09/05/2007

17:30:20

CW
3 of 7
B4

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:31:39
Diagnostic Test	Pass	17:32:07
Air Blank	RFI*	17:32:25
Air Blank	IPS**	17:32:34

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 26-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
09/05/2007

17:32:49

CM
4087
B4

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:33:39
Diagnostic Test	Pass	17:34:06
Air Blank	0.000	17:34:34
0.100 Cal Check	0.098	17:34:52
Air Blank	0.000	17:35:22
Subject Test	0.150*	17:38:30
Air Blank	0.000	17:39:02
Air Blank	0.000	17:39:29
Subject Test	XXX**	17:39:52
Air Blank	0.000	17:40:25

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000712

09/05/2007

17:40:45

CM
5077
D4

Standard Lot#: 635301I

Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:41:13
Cal Check	0.098	17:41:31
Air Blank	0.000	17:42:01



Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

CM
10-25-07
60 of 7
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 9-5-2007 TIME 1722
Intoxilyzer Serial # 80-000712 Location lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)


INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
09/05/2007 17:22:48

CM
7077
B4

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:23:16
Cal Check	0.098	17:23:35
Air Blank	0.000	17:24:04

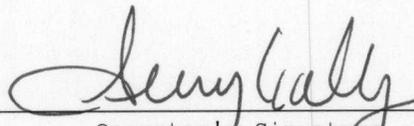

Operator's Signature

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.