



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

By  
3.17.08

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/13/08 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 2/14/08 and 3/13/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione


Sarah L. Cione

3/14/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 14, day of March, 2008.

Kathy A. Kinney  
Notary Public

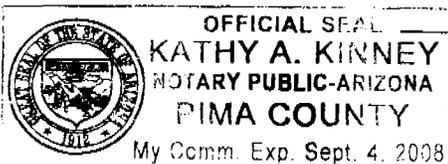

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/13/08 TIME 10:43

INTOXILYZER SERIAL # 80-000712 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
03/13/2008 10:45:05

2 of 7sc  
BY

Standard Lot#: 719802I  
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test             | g/210L | Time     |
|------------------|--------|----------|
| Air Blank        | 0.000  | 10:45:51 |
| Diagnostic Test  | Pass   | 10:46:18 |
| Air Blank        | 0.000  | 10:46:46 |
| 0.100 Cal Check  | 0.099  | 10:47:04 |
| Air Blank        | 0.000  | 10:47:34 |
| Subject Test     | 0.000  | 10:48:03 |
| Air Blank        | 0.000  | 10:48:30 |
| Five Minute Wait | ABT*   |          |
| Air Blank        | 0.000  | 10:49:22 |

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
03/13/2008 10:49:51

3 of 7sc  
154

---

Standard Lot#: 719802I  
Last Changed By: S. CLEARY #44524

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 10:50:45 |
| Diagnostic Test | Pass   | 10:51:12 |
| Air Blank       | RFI*   | 10:51:22 |
| Air Blank       | IPS**  | 10:51:32 |

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
03/13/2008 10:51:46

4 of 7sc

BA

Standard Lot#: 719802I  
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 10:52:30 |
| Diagnostic Test | Pass   | 10:52:57 |
| Air Blank       | 0.000  | 10:53:25 |
| 0.100 Cal Check | 0.099  | 10:53:43 |
| Air Blank       | 0.000  | 10:54:13 |
| Subject Test    | 0.073* | 10:57:21 |
| Air Blank       | 0.000  | 10:57:53 |
| Air Blank       | 0.000  | 10:58:20 |
| Subject Test    | XXX**  | 10:59:02 |
| Air Blank       | 0.000  | 10:59:35 |

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000712  
03/13/2008 11:00:03


5 of 75c

754

-----  
Standard Lot#: 719802I  
Last Changed By: S. CLEARY #44524  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 11:00:30 |
| Cal Check | 0.098  | 11:00:49 |
| Air Blank | 0.000  | 11:01:19 |

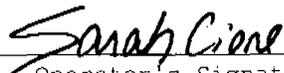  
-----  
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

104

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 3/13/08 TIME 10:43  
Intoxilyzer Serial # 80-000712 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000712  
03/13/2008 10:43:26

7 of 7sc

PS4

-----  
Standard Lot#: 719802I  
Last Changed By: S. CLEARY #44524  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 10:43:53 |
| Cal Check | 0.099  | 10:44:12 |
| Air Blank | 0.000  | 10:44:41 |

*Sarah Cione*

Operator's Signature