



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

134
6-10-2008

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/5/08 Cal 1 0.100 AC Cal 2 0.101 AC

Between the dates of 5/8/08 and 6/5/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —
The instrument was returned to service on —


Sarah Cione

Sarah L. Cione

6/9/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9, day of June, 2008.

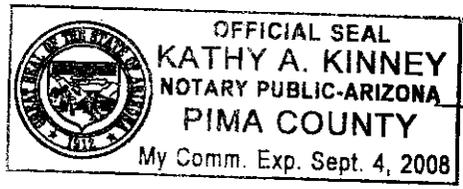

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

T54
6-10-08

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 6/5/08 TIME 9:15
INTOXILYZER SERIAL # 80-000712 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
06/05/2008 09:17:15

TS4

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:17:57
Diagnostic Test	Pass	09:18:24
Air Blank	0.000	09:18:52
0.100 Cal Check	0.100	09:19:10
Air Blank	0.000	09:19:40
Subject Test	0.000	09:20:02
Air Blank	0.000	09:20:29
Five Minute Wait	ABT*	
Air Blank	0.000	09:21:04

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
06/05/2008

09:22:03

3 of 75C
134

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:23:20
Diagnostic Test	Pass	09:23:48
Air Blank	RFI*	09:24:01
Air Blank	IPS**	09:24:16

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
06/05/2008

09:25:02

4 of 7sc

TS4

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:25:47
Diagnostic Test	Pass	09:26:14
Air Blank	0.000	09:26:41
0.100 Cal Check	0.100	09:27:00
Air Blank	0.000	09:27:29
Subject Test	0.084*	09:30:38
Air Blank	0.000	09:31:09
Air Blank	0.000	09:31:37
Subject Test	XXX**	09:31:53
Air Blank	0.000	09:32:26

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
06/05/2008 09:35:27

5 of 7 SC

B4

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:35:55
Cal Check	0.101	09:36:13
Air Blank	0.000	09:36:43

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
6-10-08

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/5/08 TIME 9:15
Intoxilyzer Serial # 80-000712 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
06/05/2008 09:15:46

7 of 75c
154

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:16:13
Cal Check	0.100	09:16:32
Air Blank	0.000	09:17:02

Sarah Cione
Operator's Signature