



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the
periodic maintenance* for the Intoxilyzer 8000 Serial Number
80-000712 maintained by the Tucson Police Department
Crime Laboratory pursuant to the requirements of the Arizona Department
of Public Safety for a quality assurance program.

By 9-2-08

A quality assurance specialist (QAS) tested the instrument on the following
dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2,
for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/2/08 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 8/1/08 and 9/2/08
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —
The instrument was returned to service on —

Sarah Cione

Sarah L. Cione


9/2/2008

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 2 day of
September, 2008.

Laura A. Santamarina

Notary Public

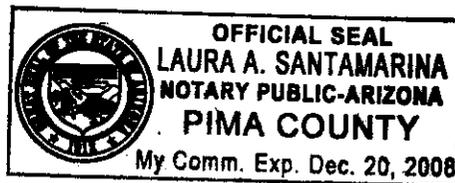

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the
information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

134
9.2.08

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/2/08 TIME 14:42

INTOXILYZER SERIAL # 80-000712 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS Calibration check greater than 31 days-

SIGNATURE Sarah Cione

2 of 7sc

ASy

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
09/02/2008

14:43:35

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:44:18
Diagnostic Test	Pass	14:44:46
Air Blank	0.000	14:45:13
0.100 Cal Check	0.099	14:45:32
Air Blank	0.000	14:46:01
Subject Test	0.000	14:46:22
Air Blank	0.000	14:46:49
Five Minute Wait	ABT*	
Air Blank	0.000	14:47:23

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7 sc

PM

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
09/02/2008 14:47:30

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	14:48:00
Air Blank	IPS**	14:48:08

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7sc

B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
09/02/2008 14:48:13

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:48:55
Diagnostic Test	Pass	14:49:23
Air Blank	0.000	14:49:50
0.100 Cal Check	0.100	14:50:09
Air Blank	0.000	14:50:38
Subject Test	0.070*	14:53:47
Air Blank	0.000	14:54:18
Air Blank	0.000	14:54:45
Subject Test	XXX**	14:55:30
Air Blank	0.000	14:56:02

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 75c

By

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
09/02/2008 14:56:25

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:56:53
Cal Check	0.098	14:57:12
Air Blank	0.000	14:57:41

Sarah Cione
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

TS4

9.2.08

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 9/2/08 TIME 14:42
Intoxilyzer Serial # 80-000712 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
09/02/2008


14:42:15

7 of 7 SC
BY

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:42:42
Cal Check	0.098	14:43:01
Air Blank	0.000	14:43:30

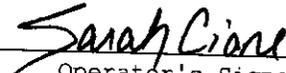

Operator's Signature