



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
1-12-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 01/07/10 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 12/10/09 and 01/07/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub
Cheryl K. Straub

01/11/10
ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 11th day of

January, 2010.

A.M. Hernandez

Notary Public



*Includes documents specified in A.A.C. R13-10-104.A.5

197
CS

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSY
1-12-2010

QA SPECIALIST Cheryl Straub AGENCY Tucson Police Department

DATE 01/07/10 TIME 7:56

INTOXILYZER SERIAL # 80-000712 LOCATION DUT. Squad / TEAM II

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

297
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000712

01/07/2010

07:59:48

Bh

Standard Lot#: 820501I

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:00:41 |
| Diagnostic Test | Pass | 08:01:09 |
| Air Blank | 0.000 | 08:01:36 |
| 0.100 Cal Check | 0.101 | 08:01:55 |
| Air Blank | 0.000 | 08:02:24 |
| Subject Test | 0.000 | 08:03:03 |
| Air Blank | ABT* | 08:03:10 |
| Air Blank | 0.000 | 08:03:38 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
01/07/2010

08:05:29

SM

Standard Lot#: 8205011
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | RFI* | 08:06:06 |
| Air Blank | IPS** | 08:06:14 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

477
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
01/07/2010

08:07:19

154

Standard Lot#: 820501I
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:08:10 |
| Diagnostic Test | Pass | 08:08:38 |
| Air Blank | 0.000 | 08:09:05 |
| 0.100 Cal Check | 0.100 | 08:09:24 |
| Air Blank | 0.000 | 08:09:53 |
| Subject Test | 0.195* | 08:13:02 |
| Air Blank | 0.000 | 08:13:34 |
| Air Blank | 0.000 | 08:14:01 |
| Subject Test | XXX** | 08:14:36 |
| Air Blank | 0.000 | 08:15:10 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

517
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
01/07/2010 08:18:08

Standard Lot#: 820501I
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

134

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 08:18:35 |
| Cal Check | 0.099 | 08:18:54 |
| Air Blank | 0.000 | 08:19:23 |



Operator's Signature

EXHIBIT G-2

697
CSTHIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4
1-12-2010

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
 DATE 01/07/10 TIME 07:56
 Intoxilyzer Serial # 80-000712 Location D4I SQUAD/TEAM II

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
 OR
 Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000712

01/07/2010

07:56:16

SH


Standard Lot#: 8205011

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 07:56:43 |
| Cal Check | 0.100 | 07:57:02 |
| Air Blank | 0.000 | 07:57:32 |

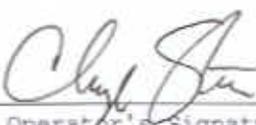

Operator's Signature