



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
6-27-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/23/11 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 5/26/11 and 6/23/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione

6/27/11


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27 day of June, 20 11.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

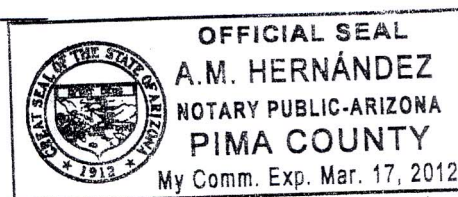

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7_{sc}

TS4
6-27-2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/23/11 TIME 926

INTOXILYZER SERIAL # 80-000712 LOCATION Team 2 - ODW

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
(☒) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
06/23/2011

09:27:47

2 of 7sc

134

Standard Lot#: 719919
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:28:47
Diagnostic Test	Pass	09:29:14
Air Blank	0.000	09:29:42
0.100 Cal Check	0.098	09:30:00
Air Blank	0.000	09:30:30
Subject Test	0.000	09:30:50
Air Blank	0.000	09:31:17
Five Minute Wait	ABT*	
Air Blank	0.000	09:31:57

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000712

06/23/2011

09:32:05

3 of 7 sc

B4

Standard Lot#: 719919

Last Changed By: S. CLEARY

#44524

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:32:53
Air Blank	IPS**	09:33:01

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000712

06/23/2011

09:33:06

Standard Lot#: 719919

Last Changed By: S. CLEARY

#44524

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:33:53
Diagnostic Test	Pass	09:34:21
Air Blank	0.000	09:34:48
0.100 Cal Check	0.098	09:35:07
Air Blank	0.000	09:35:36
Subject Test	0.120*	09:38:45
Air Blank	0.000	09:39:17
Air Blank	0.000	09:39:44
Subject Test	XXX**	09:40:05
Air Blank	0.000	09:40:38

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INFOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
06/23/2011 09:42:06

Standard Lot#: 719919
Last Changed By: S. CLEARY #44524
QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:42:33
Cal Check	0.098	09:42:52
Air Blank	0.000	09:43:22

Sarah Cione
Operator's Signature

5 of 7 sc

134

EXHIBIT G-2

6 of 7_{sc}

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
6-27-2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/23/11 TIME 926
Intoxilyzer Serial # 80-000712 Location Team 2 - ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000712

06/23/2011

09:26:08

Standard Lot#: 719919

Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:26:35
Cal Check	0.098	09:26:54
Air Blank	0.000	09:27:23

Sarah Cione
Operator's Signature