



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
6-27-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/21/12 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 5/24/12 and 6/21/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


6/22/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 22, day of June, 2012.

Ann L. Beecroft

Notary Public

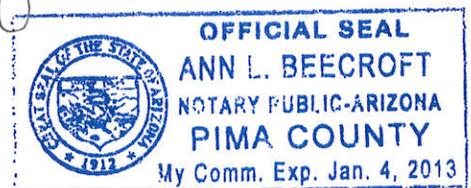

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 9<sub>sc</sub>  
TS4  
6-27-2012

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/21/12 TIME 1627

INTOXILYZER SERIAL # 80-000712 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS replaced breath hose - see attached documentation  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 9 sc  
154

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
06/21/2012 16:28:40

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:29:22
Diagnostic Test	Pass	16:29:56
Air Blank	0.000	16:30:23
0.100 Cal Check	0.098	16:30:44
Air Blank	0.000	16:31:13
Subject Test	0.000	16:31:54
Air Blank	0.000	16:32:21
Wait	ABT*	
Air Blank	0.000	16:33:24

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
06/21/2012

16:34:34

3 of 9 sc  
184

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:35:16
Diagnostic Test	Pass	16:35:50
Air Blank	0.000	16:36:17
0.100 Cal Check	0.099	16:36:38
Air Blank	0.000	16:37:08
Subject Test	RFI*	16:37:26
Air Blank	IPS**	16:37:33

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 9<sub>5c</sub>

TS4

INTOXILYZER 8000  
 Location: TUCSON PD  
 Serial Number: 80-000712  
 Core Version: 8105.48  
 06/21/2012 16:37:40

---

Standard Lot#: 23511100A2  
 Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
 TUCSON PD

---

Operator: S. CIONE #50616  
 TUCSON PD

---

Subject: SQAP, SQAP  
 DOB: 11/11/1980  
 Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:38:23
Diagnostic Test	Fail*	16:38:27
Air Blank	0.000	16:38:55

\* Voltage/Current Test

**Not a Successfully Completed Test Sequence**

5 of 9 sc

184

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
06/21/2012 16:39:10

---

Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:39:58
Diagnostic Test	Pass	16:40:32
Air Blank	0.000	16:40:59
0.100 Cal Check	0.098	16:41:20
Air Blank	0.000	16:41:50
Subject Test	0.056*	16:44:57
Air Blank	0.000	16:45:28
Air Blank	0.000	16:45:56
Subject Test	XXX**	16:46:14
Air Blank	0.000	16:46:48

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

6 of 9 sc  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
06/21/2012 16:47:01

-----  
Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	16:47:30
0.100 Cal Check	0.098	16:47:50
Air Blank	0.000	16:48:20

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7 of 9 SC

TJG  
6-27-2012

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 6/21/12 TIME 1627  
Intoxilyzer Serial # 80-000712 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

8 of 9  
5C  
734

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
06/21/2012 16:27:00

-----  
Standard Lot#: 23511100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	16:27:29
0.100 Cal Check	0.098	16:27:49
Air Blank	0.000	16:28:19

Preventative Maintenance  
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

9 of 9 SC

INTOXILYZER SERVICE RECORD

TS6  
6-27-2012

Date: 6/21/12

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-000712

Reason For Service: Fails Temperature Regulation during Diagnostics.  
Notified of issue by Officer Murphy. Breath hose not warm to  
the touch.

Date Out of Service: N/A

Calibration Check completed prior to removal from service

Date: \_\_\_\_\_

SQAP completed prior to removal from service

Date: \_\_\_\_\_

Corrective Action:  Sent to CMI, Inc.

Replaced breath hose