



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
2.19.2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/31/13 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 1/3/13 and 1/31/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


2/14/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 12/4<sup>th</sup>, day of February, 20 13.

Elma Escarcega

Notary Public

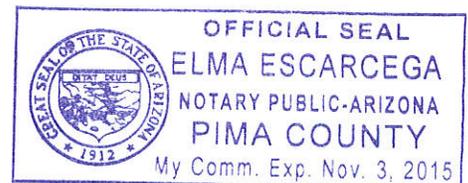

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7<sub>sc</sub>

TSG  
2.19.2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/31/13 TIME 1143

INTOXILYZER SERIAL # 80-000712 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
01/31/2013

11:47:14

2 of 7 sc  
B4

---

Standard Lot#: 09312100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:48:22
Diagnostic Test	Pass	11:48:56
Air Blank	0.000	11:49:23
0.100 Cal Check	0.099	11:49:44
Air Blank	0.000	11:50:14
Subject Test	0.000	11:51:16
Air Blank	0.000	11:51:43
Wait	ABT*	
Air Blank	0.000	11:53:04

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7<sub>sc</sub>

BY

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
01/31/2013 11:55:01

---

Standard Lot#: 09312100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	11:55:49
Air Blank	IPS**	11:55:56

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7<sub>sc</sub>

TS4

INTOXILYZER 8000  
 Location: TUCSON PD  
 Serial Number: 80-000712  
 Core Version: 8105.48  
 01/31/2013 11:59:20

Standard Lot#: 09312100A2  
 Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
 TUCSON PD

Operator: S. CIONE #50616  
 TUCSON PD

Subject: SQAP, SQAP  
 DOB: 11/11/1980  
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:00:26
Diagnostic Test	Pass	12:01:00
Air Blank	0.000	12:01:27
0.100 Cal Check	0.099	12:01:48
Air Blank	0.000	12:02:17
Subject Test	0.106*	12:05:25
Air Blank	0.000	12:05:57
Air Blank	0.000	12:06:24
Subject Test	XXX**	12:06:39
Air Blank	0.000	12:07:12

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7 sc  
TSh

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
01/31/2013 12:08:15

-----  
Standard Lot#: 09312100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:08:44
0.100 Cal Check	0.098	12:09:05
Air Blank	0.000	12:09:34

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7 sc  
TSY  
2.19.2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department

DATE 1/31/13 TIME 1143

Intoxilyzer Serial # 80-000712 Location DUIHQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0.099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 sc  
Bl

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
01/31/2013 11:43:09

-----  
Standard Lot#: 09312100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:43:38
0.100 Cal Check	0.099	11:43:58
Air Blank	0.000	11:44:28

Preventative Maintenance  
Performed -- Timer Reset