



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TSK
10-31-2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/10/13 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 9/12/13 and 10/10/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


10/21/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 21st, day of October, 2013.

Edma Escarcega

Notary Public

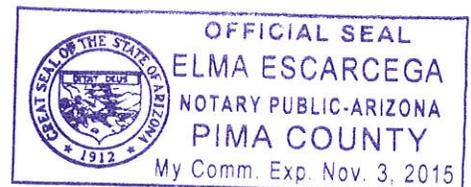

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 9 sc
TS4
10-31-2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/10/13 TIME 1056

INTOXILYZER SERIAL # 80-000712 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
10/10/2013

2 of 9_{sc}

11:04:59

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:05:46
Diagnostic Test	Pass	11:06:20
Air Blank	0.000	11:06:47
0.100 Cal Check	0.097	11:07:08
Air Blank	0.000	11:07:38
Subject Test	0.000	11:07:56
Air Blank	0.000	11:08:23
Wait	ABT*	
Air Blank	0.000	11:09:23

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
10/10/2013

11:12:16

3 of 9 sc
184

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:13:15
Diagnostic Test	Pass	11:13:49
Air Blank	0.000	11:14:17
0.100 Cal Check	0.097	11:14:38
Air Blank	0.000	11:15:07
Subject Test	RFI*	11:15:31
Air Blank	IPS**	11:15:39

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
10/10/2013

11:16:32

4 of 9 sc
Tb4

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:17:39
Diagnostic Test	Pass	11:18:13
Air Blank	0.000	11:18:40
0.100 Cal Check	0.096	11:19:01
Air Blank	0.000	11:19:31
Subject Test	0.158*	11:22:38
Air Blank	0.000	11:23:10
Air Blank	0.000	11:23:37
Subject Test	XXX**	11:24:02
Air Blank	0.000	11:24:36

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
10/10/2013 11:25:18

5 of 9 sc
TS4

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:25:47
0.100 Cal Check	0.000	11:26:00
Air Blank	0.000	11:26:37

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 9
TS4
10-31-2013

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/10/13 TIME 1056
Intoxilyzer Serial # 80-000712 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
10/10/2013 10:56:30

7 of 9
SC
Bh

Standard Lot#: 09312100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:57:44
0.100 Cal Check	0.096	10:58:05
Air Blank	0.000	10:58:35

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

8 of 9 SC
TS4
10.31.2013

Date: 10/10/13

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-000712

New Lot Number: 30412100A1

New Expiration Date: 12/1/14

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
 - Enter Target Value: 0.100
 - Enter Control Lot #: From the tank label
 - Enter Exp Date: From the tank label (mm/dd/yyyy)
 - Edit info Y/N? Y or N
 - Update Std Info? Y ENTER
 - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 096 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
10/10/2013 11:03:16

9 of 9 sc
TS4

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:03:44
0.100 Cal Check	0.096	11:04:05
Air Blank	0.000	11:04:35

Preventative Maintenance
Performed -- Timer Reset