



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
12/18/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-00072 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-11-2014 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 11-13-2014 and 12-11-2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.


Terry S. Gallegos

Terry S. Gallegos

12-17-2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 17<sup>th</sup>, day of December, 2014.

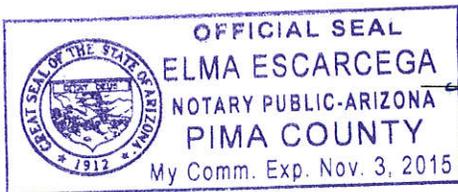

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

10 of 7  
TS  
SC  
12/18/14

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry hallegos 28847 AGENCY Tucson Police Department

DATE 12-11-2014 TIME 1533

INTOXILYZER SERIAL # 80-000712 LOCATION ODM Sgt's Office

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7  
TS6  
SC

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
12/11/2014 15:36:25

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:37:28
Diagnostic Test	Pass	15:38:05
Air Blank	0.000	15:38:32
0.100 Cal Check	0.098	15:38:53
Air Blank	0.000	15:39:22
Subject Test	0.000	15:39:36
Air Blank	0.000	15:40:03
Wait	ABT*	
Air Blank	0.000	15:40:37

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
12/11/2014 15:41:17

3of7  
TSh  
SL

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:42:14
Diagnostic Test	Pass	15:42:51
Air Blank	0.000	15:43:18
0.100 Cal Check	RFI*	15:43:32
Air Blank	IPS**	15:43:44

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 81-000712  
Core Version: 8105.48  
12/11/2014 15:44:29

4077  
TS4  
5c

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

OAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST,SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:45:17
Diagnostic Test	Pass	15:45:53
Air Blank	0.000	15:46:20
0.100 Cal Check	0.099	15:46:41
Air Blank	0.000	15:47:11
Subject Test	0.102*	15:50:18
Air Blank	0.000	15:50:50
Air Blank	0.000	15:51:18
Subject Test	XXX**	15:51:42
Air Blank	0.000	15:52:16

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
12/11/2014 15:52:37

5087  
TS4

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

5c

OAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:53:06
0.100 Cal Check	0.097	15:53:27
Air Blank	0.000	15:53:57

Preventative Maintenance  
Performed -- Timer Reset

1

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6077  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
12/18/14

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 12-11-2014 TIME 1533  
Intoxilyzer Serial # 80-000712 Location ODM Sgt's Office

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

7077  
TS4  
SC

INTOXILYZER 8000  
Instrument Initialization  
14:56:23 12/11/2014

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
12/11/2014 15:33:41

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

OAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:34:10
0.100 Cal Check	0.098	15:34:31
Air Blank	0.000	15:35:00

Preventative Maintenance  
Performed -- Timer Reset