



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
2.10.2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/5/2015 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 1/7/2015 and 2/5/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis


2/5/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 5th, day of February, 2015.

Elma Escarcega

Notary Public

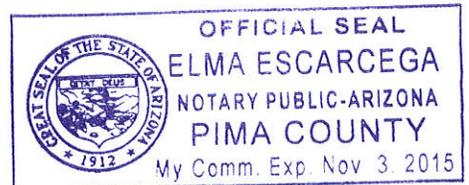

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 8
WV
TS4
2.10.2015

QA SPECIALIST WILL LOOMIS 53545 AGENCY Tucson Police Department

DATE 2/5/2015 TIME 08:57

INTOXILYZER SERIAL # 80-000712 LOCATION ODM - SGT'S OFFICE

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0.000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS RFI DETECTED DURING DIAGNOSTICS PRIOR TO
ALCOHOL-FREE SUBJECT TEST

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 800
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
02/05/2015 09:00:51

2 of 8
WLL
TBY

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

OAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Depriuation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:01:46
Diagnostic Test	RFI*	09:01:56
Air Blank	0.000	09:02:24

* RFI Detect

Not a Successfully
Completed
Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
02/05/2015 09:03:00

3 of 8
WVW

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

TS4

OAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:03:58
Diagnostic Test	Pass	09:04:35
Air Blank	0.000	09:05:02
0.100 Cal Check	0.098	09:05:23
Air Blank	0.000	09:05:53
Subject Test	0.000	09:06:09
Air Blank	0.000	09:06:36
Wait	ABT*	
Air Blank	0.000	09:07:12

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.4B
02/05/2015 09:07:44

4 of 8
WVW

TJY

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

OAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:08:33
Diagnostic Test	Pass	09:09:09
Air Blank	0.000	09:09:37
0.100 Cal Check	RFI*	09:09:52
Air Blank	IPS**	09:10:05

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
02/05/2015

09:10:45

5 of 8
WLV
Bh

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

OAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:11:41
Diagnostic Test	Pass	09:12:18
Air Blank	0.000	09:12:45
0.100 Cal Check	0.097	09:13:06
Air Blank	0.000	09:13:35
Subject Test	0.285*	09:16:43
Air Blank	0.000	09:17:16
Air Blank	0.000	09:17:43
Subject Test	XXX**	09:18:02
Air Blank	0.000	09:18:36

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

6 of 8
WV
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
02/05/2015 09:21:04

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:21:36
0.100 Cal Check	0.096	09:21:56
Air Blank	0.000	09:22:26

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7 of 8
WLL
BY
2-10-2015

QA Specialist WILL LOOMIS 53545 AGENCY Tucson Police Department
DATE 2/5/2015 TIME 08:57
Intoxilyzer Serial # 80-000712 Location ODM - SGT'S OFFICE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

8 of 8
TJ

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
02/05/2015
08:57:48

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:58:18
0.100 Cal Check	0.097	08:58:39
Air Blank	0.000	08:59:08

Preventative Maintenance
Performed -- Timer Reset