



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
8.25.2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/20/2015 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 7/23/2015 and 8/20/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

[Handwritten Signature]

Amy L. Autret


8/24/2015

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 24<sup>th</sup> day of August, 2015.

[Handwritten Signature]

Notary Public

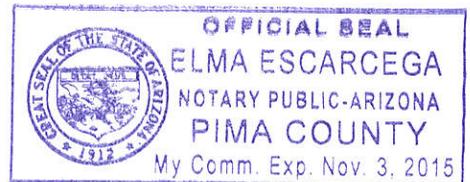

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1028  
dls  
B6  
8.25.2015

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 8/20/2015 TIME 1413

INTOXILYZER SERIAL # 80-000712 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 095 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE



20F8  
Ad  
TJ4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
08/20/2015 14:19:33

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

DAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Subject: SOAP, SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	14:20:25
Diagnostic Test	Pass	14:21:01
Air Blank	0.000	14:21:28
0.100 Cal Check	0.096	14:21:49
Air Blank	0.000	14:22:19
Subject Test	0.000	14:22:40
Air Blank	0.000	14:23:07
Wait	ABT*	
Air Blank	0.000	14:23:43

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

3 of 8  
AD  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
08/20/2015 14:25:01

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

DAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	14:25:59
Diagnostic Test	Pass	14:26:36
Air Blank	0.000	14:27:03
0.100 Cal Check	RFI*	14:27:20
Air Blank	IPS**	14:27:28

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

4 of 8  
JLB  
TSH

INTOX:LYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
08/20/2015 14:28:29

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

OAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Subject: SOAP, SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:29:26
Diagnostic Test	Pass	14:30:03
Air Blank	0.000	14:30:30
0.100 Cal Check	0.096	14:30:51
Air Blank	0.000	14:31:20
Subject Test	XXX*	14:32:33
Air Blank	0.000	14:33:07

\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

Forgot to allow  
Deficient sample to  
Time out. Repeated  
Sequence of Deficient  
Sample followed by  
Invalid sample.  
JLB

5048  
AD  
TSH

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
08/20/2015 14:34:20

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Subject: SOAP, SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	14:35:23
Diagnostic Test	Pass	14:36:00
Air Blank	0.000	14:36:27
0.100 Cal Check	0.096	14:36:48
Air Blank	0.000	14:37:17
Subject Test	0.304*	14:40:25
Air Blank	0.000	14:40:57
Air Blank	0.000	14:41:25
Subject Test	XXX**	14:41:52
Air Blank	0.000	14:42:26

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

6028  
dd  
T84

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
08/20/2015 14:44:21

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

-----  
QAS: A. AUTRET #53546  
TUCSON PD

-----  
Operator: A. AUTRET #53546  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:44:50
0.100 Cal Check	0.095	14:45:11
Air Blank	0.000	14:45:41

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

2018  
JLB  
TS4  
8.25.2015

QA Specialist Amy Antret 53546 AGENCY Tucson Police Department  
DATE 8/20/2015 TIME 1413  
Intoxilyzer Serial # 80-000712 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

8068  
[Signature]  
TS4

INTOXILYZER 8000  
Instrument Initialization  
09:41:57 08/20/2015

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000712  
Core Version: 8105.48  
08/20/2015 14:13:54

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:14:28
0.100 Cal Check	0.096	14:14:48
Air Blank	0.000	14:15:18

Preventative Maintenance  
Performed -- Timer Reset