



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSy
11-3-2015

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/15/2015 Cal 1 0.095 AC Cal 2 0.095 AC

Between the dates of 9/21/2015 and 10/15/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis


10/22/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 22nd, day of October, 2015.

Edma Escarcega

Notary Public

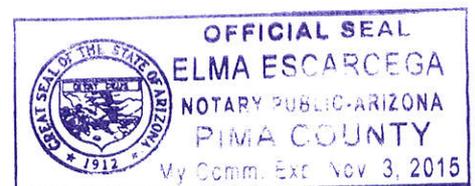

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7
WLL

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSh
11.3.2015

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 10/15/2015 TIME 16:18

INTOXILYZER SERIAL # 80-000712 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS _____

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7
wll

Bh

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
10/15/2015 16:20:58

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 16:21:46 |
| Diagnostic Test | Pass | 16:22:22 |
| Air Blank | 0.000 | 16:22:50 |
| 0.100 Cal Check | 0.095 | 16:23:10 |
| Air Blank | 0.000 | 16:23:40 |
| Subject Test | 0.000 | 16:23:59 |
| Air Blank | 0.000 | 16:24:26 |
| Wait | ABT* | |
| Air Blank | 0.000 | 16:25:01 |

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 7
WLL
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
10/15/2015 16:25:50

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

OAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 16:26:44 |
| Diagnostic Test | Pass | 16:27:21 |
| Air Blank | 0.000 | 16:27:48 |
| 0.100 Cal Check | RFI* | 16:28:00 |
| Air Blank | IPS** | 16:28:10 |

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
10/15/2015 16:29:23

4 of 7
WLL

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

TS4

OAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 16:30:50 |
| Diagnostic Test | Pass | 16:31:26 |
| Air Blank | 0.000 | 16:31:53 |
| 0.100 Cal Check | 0.095 | 16:32:14 |
| Air Blank | 0.000 | 16:32:44 |
| Subject Test | 0.219* | 16:35:51 |
| Air Blank | 0.000 | 16:36:24 |
| Air Blank | 0.000 | 16:36:51 |
| Subject Test | XXX** | 16:37:17 |
| Air Blank | 0.000 | 16:37:51 |

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 7
WLL
TBH

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
10/15/2015 16:39:01

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

OAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 16:39:30 |
| 0.100 Cal Check | 0.095 | 16:39:51 |
| Air Blank | 0.000 | 16:40:21 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7
WLL
TS4
11.3.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 10/15/2015 TIME 16:18
Intoxilyzer Serial # 80-000712 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C.
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 095 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

Feb 7 will
134

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
10/15/2015 16:18:52

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 16:19:21 |
| 0.100 Cal Check | 0.095 | 16:19:42 |
| Air Blank | 0.000 | 16:20:12 |

Preventative Maintenance
Performed - Timer Reset