



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
12-18-2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-10-15 Cal 1 0.097 AC Cal 2 0.095 AC

Between the dates of 11-12-15 and 12-10-15 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.


Amy L. Autret

12-16-15

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 16, day of December, 2015.

Notary Public

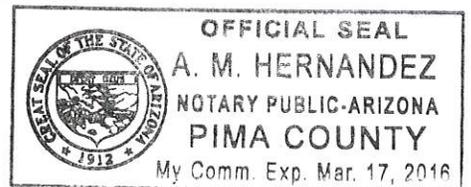

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

10F7
TS4
12.18.2015

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 12-10-15 TIME 0824

INTOXILYZER SERIAL # 80-000712 LOCATION Laboratory

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

() 1. Clock time check.

() 2. Date check.

OPERATIONAL TESTS

() 1. Alcohol-free subject test result 0. 000 AC.

() 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed


() 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

() 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

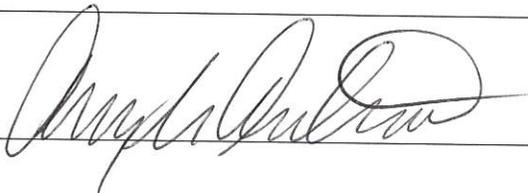

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
12/10/2015

08:27:50

2 of 7
all
TSh

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:28:49 |
| Diagnostic Test | Pass | 08:29:23 |
| Air Blank | 0.000 | 08:29:51 |
| 0.100 Cal Check | 0.096 | 08:30:11 |
| Air Blank | 0.000 | 08:30:41 |
| Subject Test | 0.000 | 08:31:02 |
| Air Blank | 0.000 | 08:31:30 |
| Wait | ABT^ | |
| Air Blank | 0.000 | 08:32:12 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
12/10/2015

08:32:45

3 of 7
all
TBH

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:34:09 |
| Diagnostic Test | Pass | 08:34:43 |
| Air Blank | 0.000 | 08:35:10 |
| 0.100 Cal Check | RFI* | 08:35:24 |
| Air Blank | IPS** | 08:35:32 |

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
12/10/2015

08:36:02

4 of 7
old
184

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:37:19 |
| Diagnostic Test | Pass | 08:37:53 |
| Air Blank | 0.000 | 08:38:20 |
| 0.100 Cal Check | 0.097 | 08:38:41 |
| Air Blank | 0.000 | 08:39:10 |
| Subject Test | 0.102* | 08:42:18 |
| Air Blank | 0.000 | 08:42:50 |
| Air Blank | 0.000 | 08:43:18 |
| Subject Test | XXX** | 08:43:38 |
| Air Blank | 0.000 | 08:44:13 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
12/10/2015 08:44:51

5 of 7
old
TBH

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:45:21 |
| 0.100 Cal Check | 0.095 | 08:45:42 |
| Air Blank | 0.000 | 08:46:11 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY


STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
T84
12-18-2015

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department
DATE 12-10-15 TIME 0824
Intoxilyzer Serial # 80-000712 Location Laboratory

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

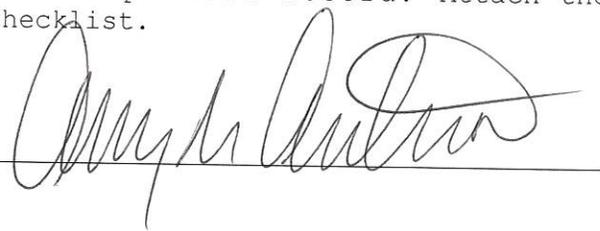SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
12/10/2015 08:24:18

7047
dd
T54

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:24:51 |
| 0.100 Cal Check | 0.097 | 08:25:12 |
| Air Blank | 0.000 | 08:25:41 |

Preventative Maintenance
Performed -- Timer Reset