



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
6.6.2016

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/26/2016 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 4/29/2016 and 5/26/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis


6/1/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 1, day of June, 2016.

A.M. Hernandez

Notary Public

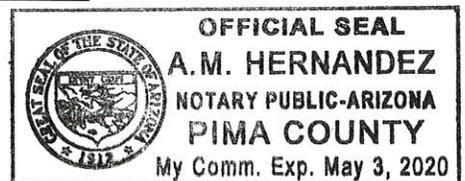

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
WLL
TSh 6-6-2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 5/26/2016 TIME 13:02

INTOXILYZER SERIAL # 80-000712 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Willie Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
05/26/2016

13:04:11

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:05:07
Diagnostic Test	Pass	13:05:41
Air Blank	0.000	13:06:08
0.100 Cal Check	0.098	13:06:29
Air Blank	0.000	13:06:58
Subject Test	0.000	13:07:12
Air Blank	0.000	13:07:40
Wait	ABT*	
Air Blank	0.000	13:08:19

*Sequence Aborted

2 of 7
WLL
TSh

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
05/26/2016

13:08:47

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:10:18
Diagnostic Test	Pass	13:10:52
Air Blank	0.000	13:11:20
0.100 Cal Check	RFI*	13:11:32
Air Blank	IPS**	13:11:39

*RFI Detect

**Improper Sample

3 of 7
WLL
TSH

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
05/26/2016

13:12:35

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:13:30
Diagnostic Test	Pass	13:14:04
Air Blank	0.000	13:14:31
0.100 Cal Check	0.097	13:14:52
Air Blank	0.000	13:15:22
Subject Test	0.074*	13:18:29
Air Blank	0.000	13:19:01
Air Blank	0.000	13:19:28
Subject Test	XXX**	13:19:39
Air Blank	0.000	13:20:13

*Deficient Sample

**Invalid Sample

Begin new deprivation period

4 of 7
WLL
TSh

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
05/26/2016 13:20:56

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:21:26
0.100 Cal Check	0.097	13:21:47
Air Blank	0.000	13:22:16

Preventative Maintenance
Performed -- Timer Reset

5 of 7
6 WLL
TBH

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
WLL
TS4 6.6.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 5/26/2016 TIME 13:02
Intoxilyzer Serial # 80-000712 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Willie Loomis

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
05/26/2016 13:02:35

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:03:04
0.100 Cal Check	0.097	13:03:25
Air Blank	0.000	13:03:55

Preventative Maintenance
Performed -- Timer Reset

7 of 7
0 WLL
TS4